

 **BD Rhapsody™ System**

Mouse TCR/BCR Next, mRNA Whole  
Transcriptome Analysis (WTA),  
BD® AbSeq, and Sample Tag  
Library Preparation Protocol

## Copyrights

No part of this publication may be reproduced, transmitted, transcribed, stored in retrieval systems, or translated into any language or computer language, in any form or by any means: electronic, mechanical, magnetic, optical, chemical, manual, or otherwise, without prior written permission from BD.

The information in this guide is subject to change without notice. BD reserves the right to change its products and services at any time. Although this guide has been prepared with every precaution to ensure accuracy, BD assumes no liability for any errors or omissions, nor for any damages resulting from the application or use of this information. BD welcomes customer input on corrections and suggestions for improvement.

## Patents and Trademarks

For US patents that may apply, see [bd.com/patents](https://bd.com/patents).

BD, the BD Logo and BD Rhapsody are trademarks of Becton, Dickinson and Company or its affiliates. All other trademarks are the property of their respective owners. © 2025 BD. All rights reserved.

## Regulatory information

For Research Use Only. Not for use in diagnostic or therapeutic procedures.

## History

Revision	Date	Change made
23-24519(01)	2024-01	Initial release.
23-24519(02)	2025-05	Updated temperatures to PCR2 conditions for TCR and BCR table and removed beads from PCR1 conditions for Targeted mRNA.

# Contents

---

<b>Introduction</b>	<b>5</b>
Protocol kits .....	5
Symbols .....	7
<b>Workflows</b>	<b>7</b>
Reverse transcription, template switching, and Exonuclease I treatment workflow .....	7
WTA library amplification workflow .....	9
BD® AbSeq/Sample Tag library amplification workflow .....	9
TCR/BCR library amplification workflow .....	10
<b>Required and recommended materials</b>	<b>11</b>
Required reagents .....	11
Recommended consumables .....	12
Equipment .....	13
Best practices .....	13
Additional documentation .....	14
Safety information .....	15
Time considerations .....	16
<b>Procedure</b>	<b>17</b>
<b>1. Reverse transcription, template switching, and Exonuclease I treatment</b>	<b>17</b>
1.1 cDNA synthesis and template switching .....	17
1.2 Denaturation and self-hybridization .....	19
1.3 TCR/BCR extension .....	20
1.4 Exonuclease I treatment .....	21
<b>2. WTA library amplification</b>	<b>22</b>
Before you begin .....	22
2.1 Random priming and extension (RPE) .....	22

2.2 RPE cleanup .....	24
2.3 RPE PCR .....	24
2.4 RPE PCR cleanup and quality check .....	26
2.5 WTA index PCR .....	27
2.6 WTA index PCR cleanup and quality check .....	29
<b>3. BD® AbSeq/Sample Tag library amplification</b> .....	<b>30</b>
3.1 BD® AbSeq/Sample Tag PCR1 .....	30
3.2 BD® AbSeq/Sample Tag PCR1 cleanup .....	31
3.3 BD® AbSeq/Sample Tag PCR1 quality check .....	32
3.4 Sample Tag PCR2 .....	33
3.5 Sample Tag PCR2 cleanup and quantify .....	34
3.6 BD® AbSeq/Sample Tag index PCR .....	35
3.7 BD® AbSeq/Sample Tag index PCR cleanup and quality check .....	37
<b>4. TCR/BCR library amplification</b> .....	<b>39</b>
Before you begin .....	39
4.1 TCR/BCR PCR1 .....	39
4.2 TCR/BCR PCR1 cleanup .....	41
4.3 TCR/BCR PCR2 .....	41
4.4 TCR/BCR PCR2 cleanup .....	43
4.5 TCR/BCR RPE .....	44
4.6 TCR/BCR RPE cleanup .....	45
4.7 TCR/BCR index PCR .....	46
4.8 TCR/BCR index PCR cleanup and quality check .....	47
<b>Sequencing</b> .....	<b>49</b>
Read requirements for libraries .....	49
Pooling libraries for sequencing .....	50
Sequencing flow cell loading and PhiX concentrations .....	53
Sequencing analysis pipeline .....	54
<b>Appendix</b> .....	<b>55</b>
Illumina Index 1 (i7) sequences .....	55
Illumina Index 2 (i5) sequence .....	55
Mouse T cell PCR1 primers .....	55
Mouse T cell PCR2 primers .....	55
Mouse B cell PCR1 primers .....	56
Mouse B cell PCR2 primers .....	56
<b>Contact Information</b> .....	<b>57</b>

## Introduction

This protocol enables high throughput single-cell transcriptome and protein analysis alongside TCR and BCR profiling of individual cells captured on the BD Rhapsody™ system, providing instructions for amplifying Illumina-compatible single-cell barcoded mRNA, TCR, BCR, BD® AbSeq, and Sample Tag libraries. The addition of a Sample Tag library enables demultiplexing of samples that were pooled before capture on the BD Rhapsody™ system.

After staining cells with BD® AbSeq and Sample Tag antibodies and partitioning and lysis of cells, cDNA is encoded on BD Rhapsody™ Enhanced Cell Capture Beads using both the 3' and 5' ends of transcripts as templates. mRNA, TCR, and BCR libraries are then amplified from these on-bead cDNA libraries using a two-step nested amplification, with TCR and BCR libraries undergoing additional random priming to capture complementarity determining regions (CDR) 1, 2, and 3, as well as framework regions (FR) 1-4. BD® AbSeq and Sample Tag libraries are amplified from the supernatant that was denatured from the beads.

### Protocol kits

Before you begin, ensure that you have the correct kits for this protocol. Matching cap colors indicate you have the correct kit, along with the catalog numbers found in the [Required and recommended materials](#) section.

#### BD Rhapsody™ cDNA Kit

Cap Color	Name	Quantity
	RT Buffer	1
	RT 0.1M DTT	1
	Reverse Transcriptase	1
	dNTP	1
	RNase Inhibitor	1
	Bead RT/PCR Enhancer	1
	10X Exonuclease I Buffer	1
	Exonuclease I	1
	Nuclease-free water	2
	Bead Resuspension Buffer	1

#### BD Rhapsody™ Enhanced Cartridge Reagent Kit V3

Cap Color	Name	Quantity
	BD Rhapsody™ HT Enhanced Cell Capture Beads V3	4
	Sample Buffer	1
	Cartridge Wash Buffer 1	1
	Cartridge Wash Buffer 2	1
	Lysis Buffer	4
	Bead Wash Buffer	1
	Waste Collection Container	4
	1M DTT	1

**BD Rhapsody™ WTA Amplification Kit**

Cap Color	Name	Quantity
	WTA Extension Primers	1
	WTA Extension Buffer	1
	WTA Extension Enzyme	1
	10 mM dNTP	1
	Nuclease-free water	3
	Bead RT/PCR Enhancer	3
	WTA Amplification Primer	1
	PCR Master Mix	1
	Universal Oligo	2
	Sample Tag PCR1 Primer	1
	Sample Tag PCR2 Primer	1
	BD® AbSeq PCR1 Primer	1
	Library Reverse Primer 1-4	1 each
	Library Forward Primer	2
	Bead Resuspension Buffer	3
	Elution Buffer	2

**BD Rhapsody™ TCR/BCR Next Amplification Kit**

Cap Color	Name	Quantity
	TCR/BCR Extension Primers	1
	TCR/BCR Extension Buffer	1
	TCR/BCR Extension Enzyme	1
	10 mM dNTP	2
	Nuclease-free water	2
	Bead RT/PCR Enhancer	1
	TSO Next	1
	TCR N1 Primer - mouse	1
	TCR N2 Primer - mouse	1
	BCR N1 Primer - mouse	1
	BCR N2 Primer - mouse	1
	PCR Master Mix	1
	TCR/BCR Universal Oligo N1	1
	TCR/BCR Universal Oligo N2	1
	Elution Buffer	2
	1M MgCl <sub>2</sub>	1
	Hybridization Buffer	4

## Symbols

The following symbols are used in this guide.


Symbol	Description
	Indicates a crucial step or stopping point.

## Workflows

### Reverse transcription, template switching, and Exonuclease I treatment workflow


#### Cell lysis in

**microwell:** The mRNA, BD<sup>®</sup> AbSeq, and Sample Tag oligos are captured by poly(T).

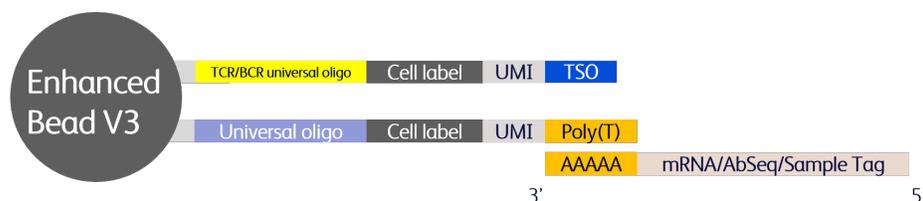


#### Reverse transcription:

Complementary DNA (cDNA) is synthesized from captured mRNA, BD<sup>®</sup> AbSeq, and Sample Tag oligos.

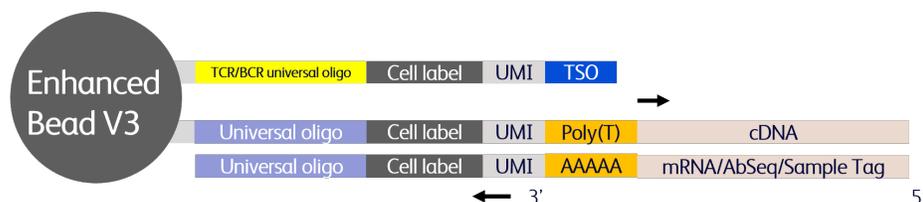

#### Template-switching

**oligo priming:** TSO is added to 3' end of cDNA on the bead.

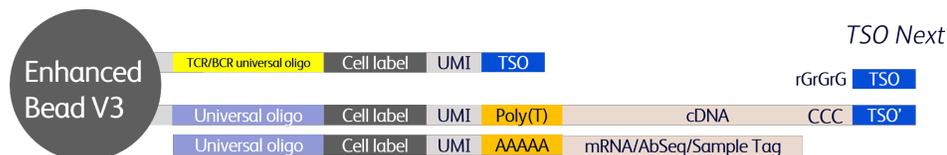

#### Denaturation


- Supernatant:** The BD<sup>®</sup> AbSeq/ Sample Tag template is denatured off the bead.




#### Important:

Save supernatant for BD<sup>®</sup> AbSeq/ Sample Tag library amplification.


- **Bead:**  
Single-stranded DNA is generated on beads to prepare for self-hybridization.

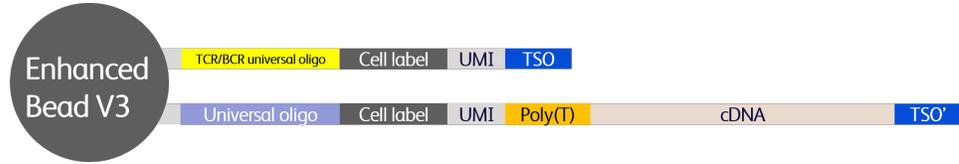

- **Self-hybridization:**  
Resuspend the beads in pre-warmed hybridization buffer, then gradually cool down to allow the TSO' to hybridize to the TSO capture strand.

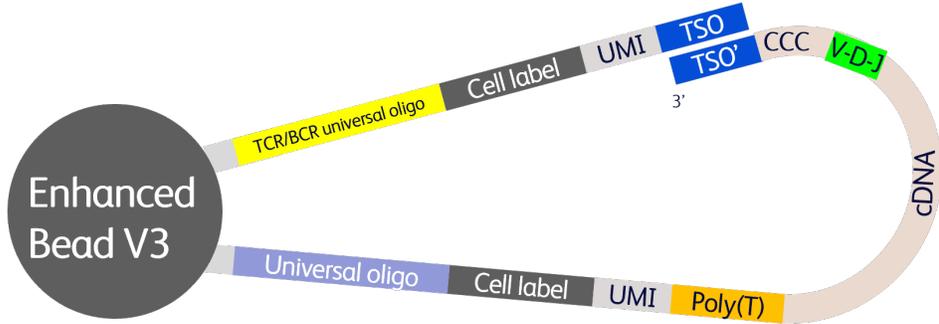


- **Extension:** Copy the TCR/BCR universal oligo onto the 3' end of cDNA strand.

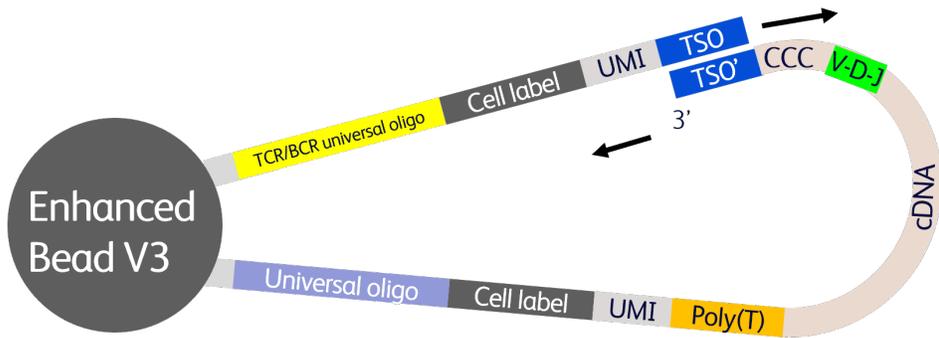

- **Exonuclease I:**  
Cleave unused oligo capture sequences off the bead.


- **Final bead layout:**  
3' sequencing template: mRNA library  
5' sequencing template: TCR/BCR library

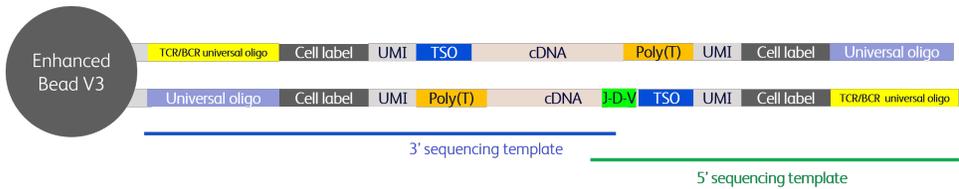

## WTA library amplification workflow


**WTA RPE:** Random priming on the bead.

Denature off the RPE product.

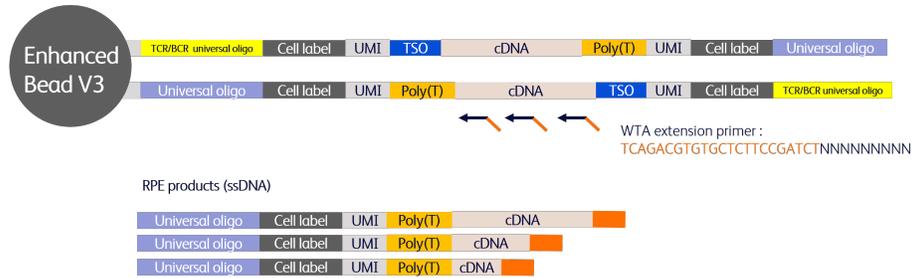


**WTA RPE**

**PCR:** Amplify the RPE product.

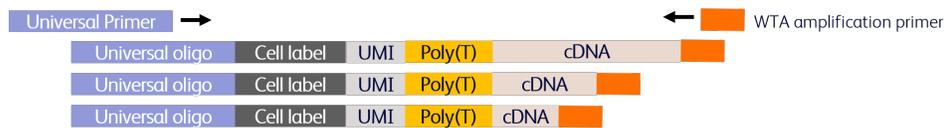


**WTA Index PCR:** Add Illumina adapters and indices.

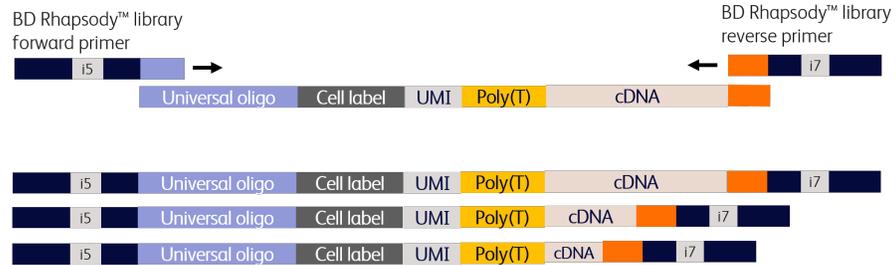

## BD® AbSeq/Sample Tag library amplification workflow


**BD® AbSeq/Sample Tag PCR1:**

BD® AbSeq/Sample Tag products are amplified together in PCR1.

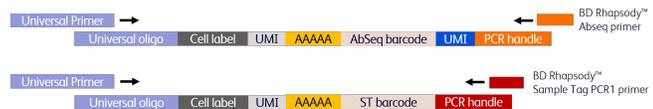


**Sample Tag**

**PCR2:** Sample Tag product is further amplified in PCR2.

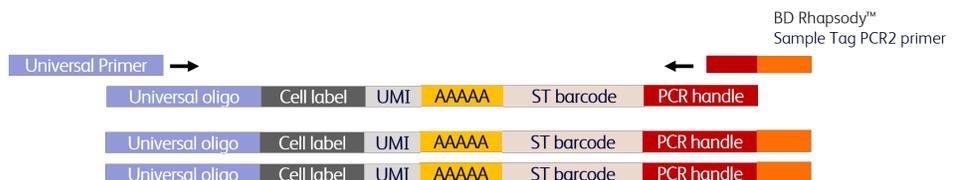

### Sample Tag Index

PCR: Add Illumina adapters and indices.

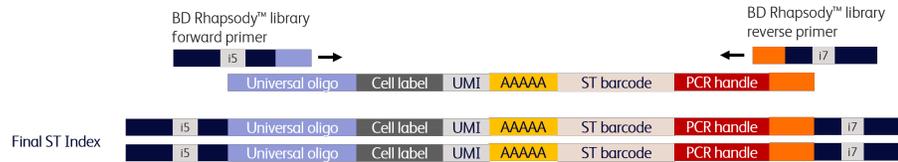

### BD® AbSeq Index

PCR: Add Illumina adapters and indices.



## TCR/BCR library amplification workflow


### TCR/BCR

PCR1: TCR/BCR universal primer and TCR/BCR N1 primer copy target region from the bead.


Amplify in solution and collect supernatant as TCR/BCR PCR1 product.

### TCR/BCR

PCR2: TCR/BCR universal N2 primer adds sequencing handle. TCR N2 or BCR N2 primer for nested PCR enrichment.

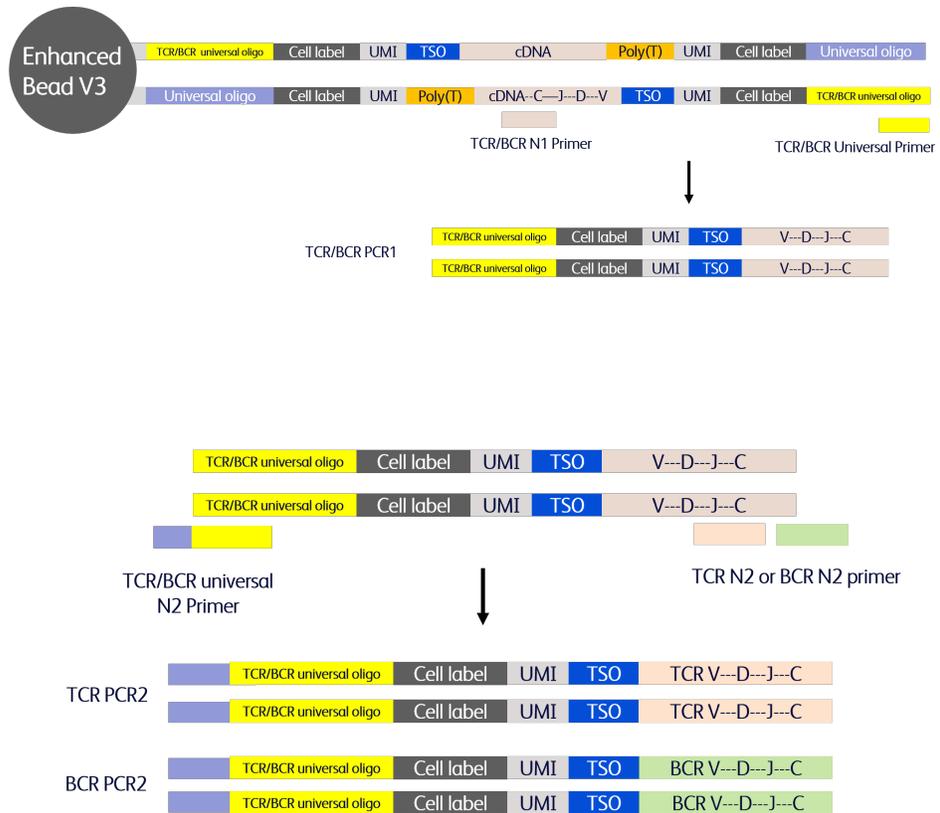


**TCR/BCR RPE for fragmentation:**  
Random priming is used for fragmentation for sequencing, full-length VDJ is assembled bioinformatically.

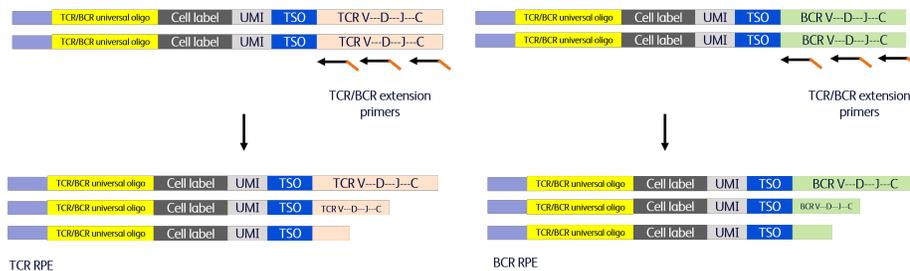

**TCR/BCR Index PCR:** Add Illumina adapters and indices.

Example: TCR indexing, BCR indexing is done separately with different indices

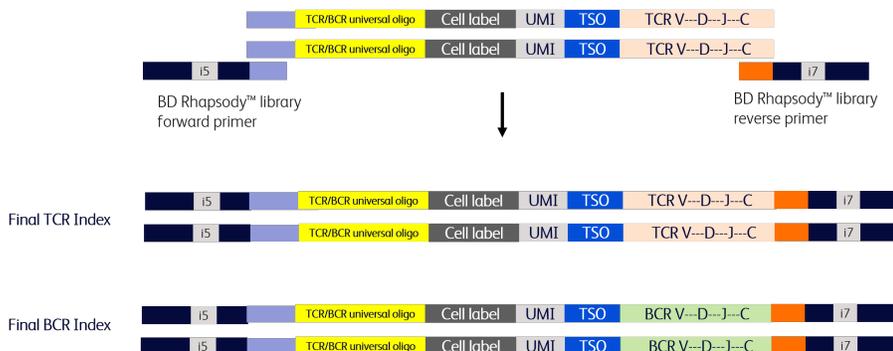

## Required and recommended materials

### Required reagents

Store the reagents at the storage temperature specified on the label.

Material	Supplier	Catalog no.
BD Rhapsody™ cDNA Kit <sup>a</sup>	BD Biosciences	633773
BD Rhapsody™ Enhanced Cartridge Reagent Kit V3 <sup>b</sup>	BD Biosciences	667052
BD Rhapsody™ WTA Amplification Kit <sup>a</sup>	BD Biosciences	633801
BD Rhapsody™ Mouse TCR/BCR Next Amplification Kit <sup>a</sup>	BD Biosciences	667059
Agencourt® AMPure® XP magnetic beads	Beckman Coulter	A63880
100% ethyl alcohol, molecular biology grade	Major supplier	–
Nuclease-free water	Major supplier	–

a. For processing more than four libraries, two orders of this catalog number are required.

b. The Enhanced Cartridge Reagent Kit V3 must be used before beginning this protocol.

**Note:** For additional indexing primers for high throughput library preparation workflows, the BD Rhapsody™ Dual Index Kit (Catalog no. 667237) is required.

## Recommended consumables

Material	Supplier	Part number/Catalog no.
Pipettes (P10, P20, P200, P1000)	Major supplier	–
Low-retention, filtered pipette tips	Major supplier	–
0.2-mL PCR 8-strip tubes	Major supplier	–
Axygen® 96-Well PCR Microplates <sup>a</sup>	Corning	PCR96HSC
Or, MicroAmp Optical 96-Well Reaction Plate <sup>a</sup>	Thermo Fisher Scientific	N8010560
MicroAmp Clear Adhesive Film	Thermo Fisher Scientific	4306311
15-mL conical tube	Major supplier	–
DNA LoBind® Tubes, 1.5 mL	Eppendorf	0030108051
DNA LoBind® Tubes, 5.0 mL	Eppendorf	0030108310
Qubit™ Assay Tubes	Thermo Fisher Scientific	Q32856
Qubit™ dsDNA HS Assay Kit	Thermo Fisher Scientific	Q32851
Agilent High Sensitivity DNA Kit	Agilent	5067-4626
Or, Agilent High Sensitivity D1000 ScreenTape	Agilent	5067-5584
Agilent High Sensitivity D1000 Reagents	Agilent	5067-5585
Or, Agilent High Sensitivity D5000 ScreenTape	Agilent	5067-5592
Agilent High Sensitivity D5000 Reagents	Agilent	5067-5593

a. Recommended for processing high throughput library preparation workflows.

## Equipment

Material	Supplier	Catalog no.
Microcentrifuge for 1.5–2.0-mL tubes	Major supplier	–
Microcentrifuge for 0.2-mL tubes	Major supplier	–
Vortexer	Major supplier	–
Digital timer	Major supplier	–
Eppendorf ThermoMixer <sup>®</sup> C	Eppendorf	538200023
6-tube magnetic separation rack for 1.5-mL tubes Or, 12-tube magnetic separation rack <sup>a</sup> Or, Invitrogen™ DynaMag™-2 magnet <sup>a</sup>	New England Biolabs New England Biolabs Thermo Fisher Scientific	S1506S S1509S 12321D
Low-profile magnetic separation stand for 0.2 mL, 8-strip tubes	V&P Scientific, Inc.	VP772F4-1
Magnetic Stand–96 <sup>b</sup>	Thermo Fisher Scientific	AM10027
Qubit™ 3.0 Fluorometer	Thermo Fisher Scientific	Q33216
Agilent <sup>®</sup> 2100 Bioanalyzer Or, Agilent <sup>®</sup> 4200 TapeStation System	Agilent Technologies Agilent Technologies	G2940CAG G2991AA
Heat block	Major supplier	–

a. Recommended for processing greater than six samples.

b. Recommended for processing high throughput library preparation workflows.

## Best practices

**⚠ Important:** The BD Rhapsody™ Enhanced Cartridge Reagent Kit V3 (Catalog no. 667052) must be used for this protocol. The BD Rhapsody™ TCR/BCR Next Amplification Kit (Catalog no. 667058) is not compatible with the BD Rhapsody™ Enhanced Cartridge Reagent Kit (Catalog no. 664887).

### Cell capture

- Ensure that the intended total cell load is 7,500–20,000. Cell loads outside this recommended range may require protocol optimization and might yield suboptimal results.
- For best results, ensure that cells have high viability before proceeding with cell capture.

### Bead handling

- When working with BD Rhapsody™ Enhanced Cell Capture Beads, use low-retention filtered tips and LoBind® Tubes.

**Never vortex the beads. Pipet-mix only.**

- Store BD Rhapsody™ Enhanced Cell Capture Beads at 4 °C. Do not freeze.
- Bring Agencourt AMPure® XP magnetic beads to room temperature (15–25 °C) before use. See the *AMPure® XP User's Guide* for information.

### Master mix preparation

- Thaw reagents (except for enzymes) at room temperature.
- Keep enzymes at –25 °C to –15 °C until ready for use.
- Return reagents to correct storage temperature as soon as possible after preparing the master mix.

### Denaturation and self-hybridization

- Remove supernatant promptly after 95 °C denaturation step (≤30 seconds after placing on magnet).
- Ensure that Hybridization Buffer is pre-heated at 80 °C for at least 20 minutes before resuspending beads in step 7 of [1.2 Denaturation and self-hybridization](#). Using cold or room temperature Hybridization Buffer may negatively impact self-hybridization efficiency.

### Supernatant handling

- Read the protocol carefully before beginning each section. Note which steps require you to keep supernatant to avoid accidentally discarding required products.
- Remove supernatants without disturbing AMPure® XP magnetic beads.
- Make and use fresh 80% ethyl alcohol within 24 hours. Adjust the volume of 80% ethyl alcohol depending on the number of libraries.

### Bead amplification

- Do not proceed to thermal cycling until each tube is gently mixed by pipette to ensure uniform bead suspension. Start the thermocycler program immediately after mixing.
- Save beads after the first amplification step ([2.1 Random priming and extension \(RPE\)](#)). They must be used again for the second bead amplification step ([4.1 TCR/BCR PCR1](#)).

### Additional documentation

- *BD Rhapsody™ Single-Cell Analysis System Instrument User Guide*
- *BD Rhapsody™ Express Single-Cell Analysis System Instrument User Guide*
- *BD Rhapsody™ System Single-Cell Labeling with BD® AbSeq Ab-Oligos for Intracellular CITE-seq Protocol*
- *BD Rhapsody™ System Single-Cell Labeling with BD Single-Cell Multiplexing Kits Protocol*

**Safety information**

For safety information, see the *BD Rhapsody™ Single-Cell Analysis Instrument User Guide* or the *BD Rhapsody™ Express Single-Cell Analysis System Instrument User Guide*.

## Time considerations

Station	Workflow		Timing	Stopping Point
Section 1: Reverse Transcription, Template Switching, and Exonuclease Treatment				
Pre-PCR Station	Step 1.1	cDNA Synthesis and Template Switching	100 min 	4 °C up to 1 week
	Step 1.2	Denature Ab-Oligos / Sample Tags	35 min	
		Self-hybridization		
	Step 1.3	TCR/BCR Extension*	30 min	
Step 1.4	Exonuclease I Treatment	50 min 	4 °C up to 3 months	
Section 2: WTA Library Amplification				
Pre-PCR Station	Step 2.1	Random Priming and Extension	210 min	
	Step 2.2	RPE Cleanup		
Post-PCR Station	Step 2.3	RPE PCR		PCR overnight
	Step 2.4	RPE PCR Cleanup and Quality Check**		<6 weeks at 4 °C or <6 months at -20 °C
	Step 2.5	WTA Index PCR	65 min	PCR overnight
	Step 2.6	WTA Index PCR Cleanup and Quality Check		<6 months at -20 °C
	(Optional) Additional WTA Index PCR Cleanup and Quality Check	25 min 	<6 months at -20 °C	
*While the TCR/BCR Extension program is running, proceed to BD® AbSeq/Sample Tag PCR1 (Step 3.1) **After RPE PCR Cleanup, proceed to TCR/BCR PCR1 (Step 4.1)				
Section 3: BD® AbSeq/Sample Tag Library Amplification				
Post-PCR Station	Step 3.1	BD® AbSeq/Sample Tag PCR1	60 min 	PCR overnight
	Steps 3.2-3.3	BD® AbSeq/Sample Tag PCR1 Cleanup and Quality Check		<24 hours at 4 °C or <6 months at -20 °C
	Step 3.4	Sample Tag PCR2	60 min 	PCR overnight
	Step 3.5	Sample Tag PCR2 Cleanup and Quantify		<24 hours at 4 °C or <6 months at -20 °C
	Step 3.6	BD® AbSeq/Sample Tag Index PCR	70 min 	PCR overnight
	Step 3.7	BD® AbSeq/Sample Tag Index PCR Cleanup and Quality Check		<6 months at -20 °C
Section 4: TCR/BCR Library Amplification				
Post-PCR Station	Step 4.1	TCR/BCR PCR1	80 min 	PCR overnight
	Step 4.2	TCR/BCR PCR1 Cleanup		<24 hours at 4 °C or <6 months at -20 °C
	Step 4.3	TCR/BCR PCR2	110 min 	PCR overnight
	Step 4.4	TCR/BCR PCR2 Cleanup		<24 hours at 4 °C or <6 months at -20 °C
	Step 4.5	TCR/BCR RPE	170 min 	PCR overnight
	Step 4.6	TCR/BCR RPE Cleanup		
	Step 4.7	TCR/BCR Index PCR		
	Step 4.8	TCR/BCR Index PCR Cleanup and Quality Check		

## Procedure

---

After staining the antibodies (as described in the *BD Rhapsody™ System Single-Cell Labeling with BD® AbSeq Ab-Oligos for Intracellular CITE-seq Protocol* or the *BD Rhapsody™ System Single-Cell Labeling with BD® Single-Cell Multiplexing Kits Protocol*), continue with this procedure.

Perform the experiment on the BD Rhapsody™ Single-Cell Analysis system following either the:

- *BD Rhapsody™ Single-Cell Analysis System Instrument User Guide*

**STOP** after the “Washing BD Rhapsody™ Enhanced Cell Capture Beads” section and follow this protocol from [Reverse transcription, template switching, and Exonuclease I treatment](#) and subsequent steps.

or

- *BD Rhapsody™ Express Single-Cell Analysis System Instrument User Guide*

**STOP** after the “Washing BD Rhapsody™ Enhanced Cell Capture Beads” section and follow this protocol from [Reverse transcription, template switching, and Exonuclease I treatment](#) and subsequent steps.

 **Important:** The BD Rhapsody™ Enhanced Cartridge Reagent Kit V3 (Catalog no. 667052) must be used for this protocol.

Ensure that the intended total cell load is between 7,500–20,000 single cells for this protocol. Cell load below or above this recommended range may not be suitable for current protocol configuration. Then proceed as described in the following procedure.

## 1. Reverse transcription, template switching, and Exonuclease I treatment

---

### 1.1 cDNA synthesis and template switching

Thaw reagents (except for the enzymes) in the BD Rhapsody™ cDNA Kit and necessary reagents from the BD Rhapsody™ TCR/BCR Next Amplification Kit at room temperature. Keep enzymes at –25 °C to –15 °C.

**Note:** This section should be performed in the pre-amplification workspace.

1. Set a thermomixer to 42 °C.
2. If performing self-hybridization on the same day, set a second thermomixer to:
  - 1,200 rpm and at 80 °C for 3 minutes.
  - 1,200 rpm and at 25 °C for 1 minute.

**Note:** The thermomixer set to 80 °C will be used as a heat block to warm the Hybridization Buffer, and then used with programmed cooling during [1.2 Denaturation and self-hybridization](#).

3. In a new 1.5-mL LoBind® tube, pipet the following reagents.

**cDNA/template switching mix**

Cap	Component	For 1 library (μL)	For 1 library with 20% overage (μL)	For 4 libraries with 20% overage (μL)	For 8 libraries with 20% overage (μL)
	RT Buffer	40.0	48.0	192.0	384.0
	dNTP	20.0	24.0	96.0	192.0
	RT 0.1 M DTT	10.0	12.0	48.0	96.0
	Bead RT/PCR Enhancer	12.0	14.4	57.6	115.2
	RNase Inhibitor	10.0	12.0	48.0	96.0
	Reverse Transcriptase	10.0	12.0	48.0	96.0
	Nuclease-free water	98.0	117.6	470.4	940.8
	<b>Total</b>	<b>200.0</b>	<b>240.0</b>	<b>960.0</b>	<b>1920.0</b>

- Gently vortex mix, briefly centrifuge, and place back on ice.
- Place the tube of washed BD Rhapsody™ Enhanced Cell Capture Beads on a 1.5-mL tube magnet for ≥2 minutes. Discard the supernatant.
- Remove the tube from the magnet and pipet 200 μL of cDNA Mix into the beads. Pipet-mix.

**Note:** Keep the prepared cDNA Mix with beads on ice until the suspension is transferred in the next step.

- Transfer the bead suspension to a new 1.5-mL LoBind® tube.
- Incubate the bead suspension on the thermomixer at 1,200 rpm and 42 °C for 30 minutes.

 **Important: Shaking is critical for this incubation.**

- While the bead suspension is still incubating at 1,200 rpm and 42 °C, pipet the following reagents in a new 1.5-mL LoBind® tube.

**Note:** Prepare the TSO Mix approximately within 2 minutes before the 30 minutes incubation at 42 °C is finished.

 **Important: Use immediately.**

**TSO Mix**

Cap	Component	For 1 library (μL)	For 1 library with 20% overage (μL)	For 4 libraries with 20% overage (μL)	For 8 libraries with 20% overage (μL)
	TSO Next	6.0	7.2	28.8	57.6
	1M MgCl <sub>2</sub>	2.0	2.4	9.6	19.2

Cap	Component	For 1 library (μL)	For 1 library with 20% overage (μL)	For 4 libraries with 20% overage (μL)	For 8 libraries with 20% overage (μL)
	<b>Total</b>	<b>8.0</b>	<b>9.6</b>	<b>38.4</b>	<b>76.8</b>

10. Gently vortex mix, briefly centrifuge, and keep on ice.
11. Add 8 μL of TSO Mix to the reaction, gently pipet-mix, and incubate on the thermomixer for another 30 minutes at 1,200 rpm and 42 °C.

 **Important:** If you are performing self-hybridization on the same day, complete steps 3 and 4 from [1.2 Denaturation and self-hybridization](#) now.

12. Place the bead suspension on the 1.5-mL tube magnet until the solution is clear ( $\leq 1$  minute). Discard the supernatant.
13. Remove the tube from the magnet and pipet 75 μL of Elution Buffer into the tube. Pipet-mix. Place on ice.

**OPTIONAL STOPPING POINT:** BD Rhapsody™ Enhanced Cell Capture Beads can be stored up to 7 days at 2–8 °C after template switching.

14. If using the *BD Rhapsody™ Single-Cell Analysis System Instrument User Guide*, view the BD Rhapsody™ Scanner image analysis to see if the analysis metrics passed.

## 1.2 Denaturation and self-hybridization

Thaw the Hybridization Buffer at room temperature.

1. Set a heat block to 95 °C.
2. Program the thermomixer.
  - a. 1,200 rpm and at 80 °C for 3 minutes.
  - b. 1,200 rpm and at 25 °C for 1 minute.

**Note:** If you performed cDNA synthesis on the same day, this is the same thermomixer from step 2 of [1.1 cDNA synthesis and template switching](#), and the thermomixer is already programmed.

 **Important:** Confirm "Time Mode" on the thermomixer is set to "Temperature Control" to ensure that the 25 °C temperature is reached before the 1 minute at 25 °C (step b) begins.

3. Prepare Hybridization Buffer for self-hybridization. Aliquot 1.2 mL Hybridization Buffer into a new 1.5-mL LoBind® tube and place the tube in the pre-heated 80 °C thermomixer (from step 2) without shaking.
4. Keep the tube of Hybridization Buffer in the 80 °C thermomixer until ready to use, at least 20 minutes before resuspending beads in step 7 of this section.
5. To denature, incubate the tube in the following order:
  - a. Ensure that the beads are resuspended. Pipet-mix to resuspend, if needed.
  - b. Incubate the tube at 95 °C in a heat block for 5 minutes.
  - c. Immediately after the completion of the 95 °C incubation, slightly open the lid of the tube to release air pressure within the tube.

6. Immediately place the tube on the magnet for  $\leq 30$  seconds until clear. **Keep the supernatant.** Remove the supernatant and transfer to a new 1.5-mL LoBind® tube. This contains the BD® AbSeq/Sample Tag supernatant products.

**Note:** To minimize BD® AbSeq/Sample Tag contamination in the TCR/BCR and WTA libraries, ensure that all liquid is removed from the tube. Keep the supernatant tube at 4 °C until ready to proceed to [3.1 BD® AbSeq/Sample Tag PCR1](#).

7. Resuspend the beads in 1.0 mL of pre-heated 80 °C Hybridization Buffer, and immediately place in the pre-programmed thermomixer from step 2. Start the program.

**Note:** Incubation will take approximately 25 minutes.

8. After the hybridization step, place tube on ice for at least 1 minute while TCR/BCR Extension Mix is being prepared.

### 1.3 TCR/BCR extension

Thaw the reagents for TCR/BCR Extension at room temperature. Keep TCR/BCR Extension enzyme at  $-25$  °C to  $-15$  °C.

1. Set a thermomixer to 37 °C.
2. Ensure all reagents other than the TCR/BCR Extension enzyme are at room temperature.
3. In a new 1.5-mL LoBind® tube, pipet the following reagents.

#### TCR/BCR Extension Mix

Cap	Component	For 1 library (µL)	For 1 library with 20% overage (µL)	For 4 libraries with 20% overage (µL)	For 8 libraries with 20% overage (µL)
○	TCR/BCR Extension Buffer	20	24	96	192
●	dNTP	20	24	96	192
○	TCR/BCR Extension Enzyme	10	12	48	96
○	Nuclease-free water	150	180	720	1440
	<b>Total</b>	<b>200</b>	<b>240</b>	<b>960</b>	<b>1920</b>

4. Gently vortex mix, briefly centrifuge, and keep at room temperature.
5. Briefly spin the tube with the bead suspension.
6. Place the tube of BD Rhapsody™ Enhanced Cell Capture Beads on a 1.5-mL tube magnet for  $\leq 2$  minutes. Discard the supernatant.
7. Remove the tubes from magnet and resuspend using 200 µL of TCR/BCR Extension Mix. Pipet-mix.
8. Incubate the bead suspension on a thermomixer at 1,200 rpm and 37 °C for 30 minutes.

**Note:** During TCR/BCR Extension incubation, begin BD® AbSeq/Sample Tag PCR1. See [3.1 BD® AbSeq/Sample Tag PCR1](#). You can leave the BD® AbSeq/Sample Tag PCR1 reaction in the thermocycler

when complete. TCR/BCR PCR1 can be performed after [2.4 RPE PCR cleanup and quality check](#). All PCR1 product purification (TCR/BCR and BD<sup>®</sup> AbSeq/Sample Tag) can be done at the same time.

- Briefly spin the tube with the beads suspension and place the tube on ice.

## 1.4 Exonuclease I treatment

Thaw reagents for Exonuclease I treatment at room temperature. Keep Exonuclease I enzyme at  $-25\text{ }^{\circ}\text{C}$  to  $-15\text{ }^{\circ}\text{C}$ .

- Set one thermomixer to  $37\text{ }^{\circ}\text{C}$  and a heat block to  $80\text{ }^{\circ}\text{C}$ .
- In a new 1.5-mL LoBind<sup>®</sup> tube, pipet the following reagents.

### Exonuclease I Mix

Cap	Component	For 1 library (μL)	For 1 library with 20% overage (μL)	For 4 libraries with 20% overage (μL)	For 8 libraries with 20% overage (μL)
	10X Exonuclease I Buffer	20	24	96	192
	Exonuclease I	10	12	48	96
	Nuclease-free water	170	204	816	1632
	<b>Total</b>	<b>200</b>	<b>240</b>	<b>960</b>	<b>1920</b>

- Gently vortex mix, briefly centrifuge, and keep at room temperature.
- Place the tube of BD Rhapsody<sup>™</sup> Enhanced Cell Capture Beads with TCR/BCR Extension Mix on a 1.5-mL tube magnet for  $\leq 1$  minute. Discard the supernatant.
- Remove the tube from the magnet and pipet 200  $\mu\text{L}$  Exonuclease I Mix into the tube. Pipet-mix.
- Incubate the bead suspension on thermomixer at 1,200 rpm and  $37\text{ }^{\circ}\text{C}$  for 30 minutes.
- Incubate the bead suspension in the heat block at  $80\text{ }^{\circ}\text{C}$  for 20 minutes.
- Place the tube on ice for  $\sim 1$  minute.
- Briefly spin the tube with the bead suspension.
- Place the tube on the magnet for  $\leq 1$  minute until clear. Discard the supernatant.
- Remove the tube from the magnet and pipet 200  $\mu\text{L}$  of cold Bead Resuspension Buffer into the tube. Pipet-mix.

 **STOPPING POINT:** Exonuclease I-treated beads can be stored at  $2\text{--}8\text{ }^{\circ}\text{C}$  for up to 3 months.

- Proceed to library preparation.

## 2. WTA library amplification

### Before you begin

- Obtain Exonuclease I-treated and inactivated BD Rhapsody™ Enhanced Cell Capture Beads.
- Thaw reagents in the BD Rhapsody™ WTA Amplification Kit at room temperature (15–25 °C), then immediately place on ice.

### 2.1 Random priming and extension (RPE)

This section describes how to generate random priming products. First, random primers are hybridized to the cDNA on the BD Rhapsody™ Enhanced Cell Capture Beads, followed by extension with an enzyme.

**Note:** Perform this procedure in the pre-amplification workspace.

1. Set a heat block to 95 °C and set two thermomixers to 37 °C and 25 °C, respectively.
2. In a new 1.5-mL LoBind® tube, pipet the following reagents.

#### Random Primer Mix

Cap	Component	For 1 library (µL)	For 1 library with 20% overage (µL)	For 4 libraries with 20% overage (µL)	For 8 libraries with 20% overage (µL)
●	WTA Extension Buffer	20.0	24.0	96.0	192.0
●	WTA Extension Primers	20.0	24.0	96.0	192.0
○	Nuclease-free water	134.0	160.8	643.2	1286.4
	<b>Total</b>	<b>174.0</b>	<b>208.8</b>	<b>835.2</b>	<b>1670.4</b>

3. Pipet-mix the Random Primer Mix and keep at room temperature.
 

**Note:** Use the entire sample of beads. Sub-sampling beads is not recommended for TCR/BCR FL + WTA + BD® AbSeq + Sample Tag combination assays.
4. Resuspend the Exonuclease I-treated BD Rhapsody™ Enhanced Cell Capture Beads with a pipette.
5. Place the tube of Exonuclease I-treated beads in Bead Resuspension Buffer on the 1.5-mL magnet for 2 minutes. Remove the supernatant.
6. Briefly centrifuge the tube, then place the tube on a 1.5-mL magnet for 2 minutes. Discard the supernatant.
7. Remove the tube with the BD Rhapsody™ Enhanced Cell Capture Beads from the magnet, and use a low-retention tip to pipet 174 µL of Random Primer Mix into the tube. Pipet-mix 10 times to resuspend the beads.
8. Incubate the tube in the following order:
  - a. 95 °C in a heat block (no shaking) for 5 minutes.
  - b. Thermomixer at 1,200 rpm and at 37 °C for 5 minutes.

- c. Thermomixer at 1,200 rpm and at 25 °C for 15 minutes.
9. Briefly centrifuge the tube and keep it at room temperature.
10. In a new 1.5-mL LoBind<sup>®</sup> tube, pipet the following reagents.

#### Extension Enzyme Mix

Cap	Component	For 1 library (μL)	For 1 library with 20% overage (μL)	For 4 libraries with 20% overage (μL)	For 8 libraries with 20% overage (μL)
	dNTP	8.0	9.6	38.4	76.8
	Bead RT/PCR Enhancer	12.0	14.4	57.6	115.2
	WTA Extension Enzyme	6.0	7.2	28.8	57.6
	<b>Total</b>	<b>26.0</b>	<b>31.2</b>	<b>124.8</b>	<b>249.6</b>

11. Pipet-mix the Extension Enzyme Mix.
12. Pipet 26 μL of the Extension Enzyme Mix into the sample tube containing the beads (for a total volume of 200 μL) and keep at room temperature until ready.
13. Program the thermomixer.
  - a. 1,200 rpm and at 25 °C for 10 minutes.
  - b. 1,200 rpm and at 37 °C for 15 minutes.
  - c. 1,200 rpm and at 45 °C for 10 minutes.
  - d. 1,200 rpm and at 55 °C for 10 minutes.

 **Important:** Confirm “Time Mode” is set to “Time Control” before the program begins.

14. Place the tube from step 12 in the thermomixer. Remove the tube after the program is finished.
15. Place the tube in a 1.5-mL tube magnet and discard the supernatant.
16. Remove the tube from the magnet and resuspend the beads in 205 μL of Elution Buffer using a P200 pipette.
17. To denature the random priming products off the beads, pipet to resuspend the beads.
  - a. Incubate the sample at 95 °C in a heat block for 5 minutes (no shaking).
  - b. Immediately after completing the 95 °C incubation, slightly open the lid of the tube to release air pressure within the tube.
 

**Note:** Do not incubate for more than 5 minutes.
18. Place the tube in a 1.5-mL tube magnet.
  - a. Immediately transfer 200 μL of the supernatant containing the Random Primer Extension Product (RPE Product) to a new 1.5-mL LoBind<sup>®</sup> tube and keep at room temperature.
  - b. Proceed to [2.2 RPE cleanup](#) in the following section.

19. Pipet 200  $\mu\text{L}$  of cold Bead Resuspension Buffer to the tube with leftover beads. Gently resuspend the beads by pipet-mixing only. Do not vortex. Store the beads on ice or at 4 °C in the pre-amplification workspace until needed.

 **Important:** These beads will be used for target specific amplification during [4. TCR/BCR library amplification](#). **DO NOT THROW AWAY!**

## 2.2 RPE cleanup

This section describes how to perform a single-sided AMPure® cleanup, which removes primer dimers and other small molecular weight byproducts. The final product is purified single-stranded DNA.

**Note:** Perform the purification in the pre-amplification workspace.

1. In a new 15-mL conical tube, prepare 10 mL of fresh 80% (v/v) ethyl alcohol by pipetting 8 mL of absolute ethyl alcohol to 2 mL of nuclease-free water. Vortex the tube for 10 seconds.  
**Note:** Make fresh 80% ethyl alcohol and use within 24 hours. The 80% ethyl alcohol volume should be adjusted depending on the number of libraries.
2. Bring Agencourt® AMPure® XP magnetic beads to room temperature. Vortex the AMPure® XP magnetic beads at high speed for 1 minute until the beads are fully resuspended.
3. Pipet 320  $\mu\text{L}$  of AMPure® XP magnetic beads into the tube containing the 200  $\mu\text{L}$  of RPE product supernatant. Pipet-mix at least 10 times, then briefly centrifuge.
4. Incubate at room temperature for 10 minutes.
5. Place the tube on the 1.5-mL tube magnet for 5 minutes. Discard the supernatant.
6. Keeping the tube on the magnet, gently add 1 mL of fresh 80% ethyl alcohol into the tube and incubate for 30 seconds. Discard the supernatant.
7. Repeat step 6 for a total of two washes.
8. Keeping the tube on the magnet, use a small-volume pipette to remove and discard any residual supernatant from the tube.
9. Air-dry the beads at room temperature for 5 minutes or until the beads no longer look glossy.
10. Remove the tube from the magnet and resuspend the bead pellet in 40  $\mu\text{L}$  of Elution Buffer. Pipet-mix the suspension at least 10 times until the beads are fully suspended.
11. Incubate the sample at room temperature for 2 minutes. Briefly centrifuge the tube to collect the contents at the bottom.
12. Place the tube on the magnet until the solution is clear, usually ~30 seconds.
13. Pipet the eluate (~40  $\mu\text{L}$ ) to a new PCR tube. This is the purified RPE product.
14. Keep on ice until ready to proceed with [2.3 RPE PCR](#).

## 2.3 RPE PCR

This section describes how to generate more RPE product through PCR amplification, so that there are multiple copies of each random-primed molecule.

1. In the pre-amplification workspace, in a new 1.5-mL LoBind<sup>®</sup> tube, pipet the following components.

#### RPE PCR Mix

Cap	Component	For 1 library (μL)	For 1 library with 20% overage (μL)	For 4 libraries with 20% overage (μL)	For 8 libraries with 20% overage (μL)
○	PCR MasterMix	60	72	288	576
○	Universal Oligo	10	12	48	96
○	WTA Amplification Primer	10	12	48	96
	<b>Total</b>	<b>80</b>	<b>96</b>	<b>384</b>	<b>768</b>

2. Add 80 μL of the RPE PCR Mix to the tube with the 40 μL of purified RPE product. Pipet-mix 10 times to create the RPE PCR Reaction Mix.
3. Split the mix into two PCR tubes with 60 μL of RPE PCR Reaction Mix per tube.
4. Bring the RPE PCR Reaction Mix to the post-amplification workspace and run the following PCR program.

#### RPE PCR program

Step	Cycles	Temperature	Time
Hot start	1	95 °C	3 min
Denaturation	Refer to the Recommended number of PCR cycles table <sup>a</sup>	95 °C	30 seconds
Annealing		60 °C	1 min
Extension		72 °C	1 min
Final extension	1	72 °C	2 min
Hold	1	4 °C	∞

a. Recommended PCR cycles might need to be optimized for different cell types and cell number.

#### Recommended number of PCR cycles

Number of cells in RPE	Recommended PCR cycles for PBMCs
7,500	13
10,000	12
20,000	11

5. When the RPE PCR reaction is complete, briefly centrifuge to collect the contents at the bottom of the tubes.

## 2.4 RPE PCR cleanup and quality check

This section describes how to perform a single-sided AMPure® cleanup to remove unwanted small molecular weight products from the RPE products. The final product is purified double-stranded DNA (~200–2,000 bp).

**Note:** Perform the purification in the post-amplification workspace.

1. Combine the two 60- $\mu$ L RPE PCR reactions into a new 1.5-mL tube.
2. Briefly centrifuge the tube with the RPE PCR product.

 **Important:** It is critical for the final volume to be exactly 120  $\mu$ L to achieve the appropriate size selection of the purified RPE PCR product.

3. In a new 15-mL conical tube, prepare 5 mL of fresh 80% (v/v) ethyl alcohol by pipetting 4 mL of absolute ethyl alcohol to 1 mL of nuclease-free water. Vortex the tube for 10 seconds.

**Note:** Make fresh 80% ethyl alcohol and use within 24 hours. The 80% ethyl alcohol volume should be adjusted depending on the number of libraries.


4. Bring AMPure® XP magnetic beads to room temperature. Vortex the AMPure® XP magnetic beads at high speed for 1 minute until the beads are fully resuspended.
5. Pipet 96  $\mu$ L of AMPure® XP magnetic beads into the tube containing 120  $\mu$ L of RPE PCR product. Pipet-mix at least 10 times, then briefly centrifuge the samples.
6. Incubate at room temperature for 5 minutes.
7. Place the 1.5-mL LoBind® tube on the magnet for 5 minutes. Discard the supernatant.
8. Keeping the tube on the magnet, gently add 200  $\mu$ L of fresh 80% ethyl alcohol into the tube and incubate for 30 seconds. Discard the supernatant without disturbing the beads.
9. Repeat step 8 once for a total of two washes.
10. Keeping the tube on the magnet, use a small-volume pipette to remove any residual supernatant from the tube.
11. Air-dry the beads at room temperature for 3 minutes.
12. Remove the tube from the magnet and pipet 40  $\mu$ L of Elution Buffer into the tube. Pipet-mix the suspension at least 10 times until the beads are fully suspended.
13. Incubate the sample at room temperature for 2 minutes. Briefly centrifuge the tube to collect the contents at the bottom.
14. Place the tube on the magnet until the solution is clear, usually ~30 seconds.
15. Pipet the eluate (~40  $\mu$ L) into a new 1.5-mL LoBind® tube. The RPE PCR product is ready for WTA Index PCR.

 **STOPPING POINT:** The RPE PCR libraries can be stored at –20 °C for up to 6 months or 4 °C for up to 6 weeks.

16. Perform quality control of the RPE PCR products by quantifying 2  $\mu$ L of the final sequencing library with a Qubit™ Fluorometer using the Qubit™ dsDNA HS Assay and with the Agilent 2100 Bioanalyzer using the Agilent High Sensitivity DNA Kit.

- The expected concentration from the Qubit™ Fluorometer is ~0.5 to 10 ng/μL.
- The Bioanalyzer trace should show a broad peak from ~150 to 2,000 bp. Use the concentration from 200 to 600 bp to calculate how much template to add into Index PCR. Refer to the blue-boxed regions in the sample trace images in Figure 1.

**Figure 1** Sample Bioanalyzer High Sensitivity DNA trace - RPE PCR product trace

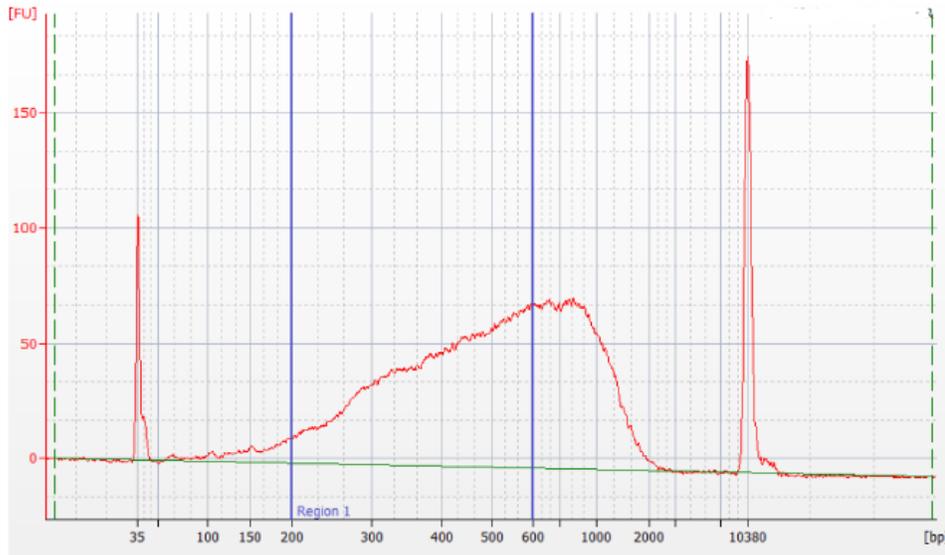

## 2.5 WTA index PCR

This section describes how to generate mRNA libraries compatible with the Illumina sequencing platform, by adding full-length Illumina sequencing adapters and indices through PCR.

**Note:** Perform this procedure in the post-amplification workspace.

- Dilute the RPE PCR products from [2.4 RPE PCR cleanup and quality check](#) with nuclease-free water such that the concentration of the 200–600 bp peak is 2 nM. If the product concentration is <2 nM, do not dilute. Continue to the next step.

**Example:** If the Bioanalyzer measurement of the 200–600 bp peak is 6 nM, then dilute the sample threefold with nuclease-free water to 2 nM.

- In a new 1.5-mL tube, pipet the following components.

**WTA Index PCR Mix**

Cap	Component	For 1 library (μL)	For 1 library with 20% overage (μL)	For 4 libraries with 20% overage (μL)	For 8 libraries with 20% overage (μL)
○	PCR MasterMix	25	30	120	240
●	Library Forward Primer	5	6	24	48
●	Library Reverse Primer 1–4 <sup>a</sup>	5	6	–	–
○	Nuclease-free water	5	6	24	48
	<b>Total</b>	<b>40</b>	<b>48</b>	<b>168</b>	<b>336</b>

a. For more than one library, use different Library Reverse Primers for each library.

- Gently vortex mix, briefly centrifuge, and place back on ice.
- In a new 0.2-mL PCR tube, combine WTA Index PCR Mix with diluted RPE PCR products as follows:
  - For one sample, combine 40 μL of WTA Index PCR Mix with 10 μL of 2 nm of RPE PCR product.
  - For multiple samples, combine 35 μL of WTA Index PCR Mix with 5 μL of Library Reverse Primer and 10 μL of 2 nm of RPE PCR products.
- Pipet-mix 10 times.
- Run the following PCR program.

**Index PCR conditions for WTA**

Step	Cycles	Temperature	Time
Hot start	1	95 °C	3 min
Denaturation	Refer to the Recommended number of PCR cycles table <sup>a</sup>	95 °C	30 seconds
Annealing		60 °C	30 seconds
Extension		72 °C	30 seconds
Final extension	1	72 °C	1 min
Hold	1	4 °C	∞

a. Cycle number varies based on the concentration of the RPE PCR product.

**Recommended number of PCR cycles**

Concentration of diluted RPE PCR products	Recommended number of PCR cycles
1 to <2 nm	9
2 nm	8

If the concentrations of diluted RPE PCR products are <1 nm, additional PCR cycles might be needed.

 **STOPPING POINT:** The PCR can run overnight.

- When the WTA Index PCR is complete, briefly centrifuge to collect the contents at the bottom of the tubes.

## 2.6 WTA index PCR cleanup and quality check

This section describes how to perform a single-sided AMPure<sup>®</sup> cleanup for Illumina sequencing. The final product is purified double-stranded DNA with full-length Illumina adapter sequences.

**Note:** Perform the purification in the post-amplification workspace.

- Add 60  $\mu\text{L}$  of nuclease-free water to the WTA Index PCR product for a final volume of 110  $\mu\text{L}$ .
- Transfer 100  $\mu\text{L}$  of WTA Index PCR product into a new 0.2-mL PCR tube.
- In a new 5.0-mL LoBind<sup>®</sup> tube, prepare 5 mL fresh 80% (v/v) ethyl alcohol by combining 4 mL absolute ethyl alcohol, molecular biology grade, with 1 mL of nuclease-free water. Vortex the tube for 10 seconds to mix.
 


**Note:** Make fresh 80% ethyl alcohol, and use it within 24 hours.
- Bring AMPure<sup>®</sup> XP beads to room temperature and vortex at high speed for 1 minute until the beads are fully resuspended.
- Add 65  $\mu\text{L}$  of AMPure<sup>®</sup> XP magnetic beads to the 0.2-mL PCR tube from step 2.
- Pipet-mix 10 times and incubate at room temperature for 5 minutes.
- Place the tube on the strip tube magnet for 3 minutes. Discard the supernatant.
- Keeping the tube on the magnet, gently add 200  $\mu\text{L}$  of fresh 80% ethyl alcohol into the tube and incubate for 30 seconds. Discard the supernatant.
- Repeat step 8 for a total of two washes.
- Keeping the tubes on the magnet, use a small-volume pipette to remove any residual supernatant from the tube.
- Leave the tubes open on the magnet to dry the AMPure<sup>®</sup> XP magnetic beads at room temperature for ~1 minute. Do not over-dry the AMPure<sup>®</sup> XP magnetic beads.
- Pipet 30  $\mu\text{L}$  of Elution Buffer into the tubes and pipet-mix to completely resuspend the AMPure<sup>®</sup> XP magnetic beads.
- Incubate the samples at room temperature for 2 minutes.
- Briefly centrifuge the tubes to collect the contents at the bottom.
- Place the tubes on the magnet until the solution is clear, usually ~30 seconds.
- Pipet the eluate (~30  $\mu\text{L}$ ) into new 1.5-mL LoBind<sup>®</sup> tubes. The WTA Index PCR eluate is the final sequencing libraries.

 **STOPPING POINT:** The Index PCR libraries can be stored at  $-20\text{ }^{\circ}\text{C}$  for up to 6 months until sequencing.

- Perform quality control of the Index PCR libraries by quantifying 2  $\mu\text{L}$  of the final sequencing library with a Qubit<sup>™</sup> Fluorometer using the Qubit<sup>™</sup> dsDNA HS Assay and with the Agilent 2100 Bioanalyzer using the Agilent High Sensitivity DNA Kit.

- a. The expected concentration from the Qubit™ Fluorometer is >1 ng/μL.
- b. The Bioanalyzer trace should show a peak from ~300 to 2,000 bp.

**Figure 2** Sample Bioanalyzer High Sensitivity DNA trace - WTA index PCR product

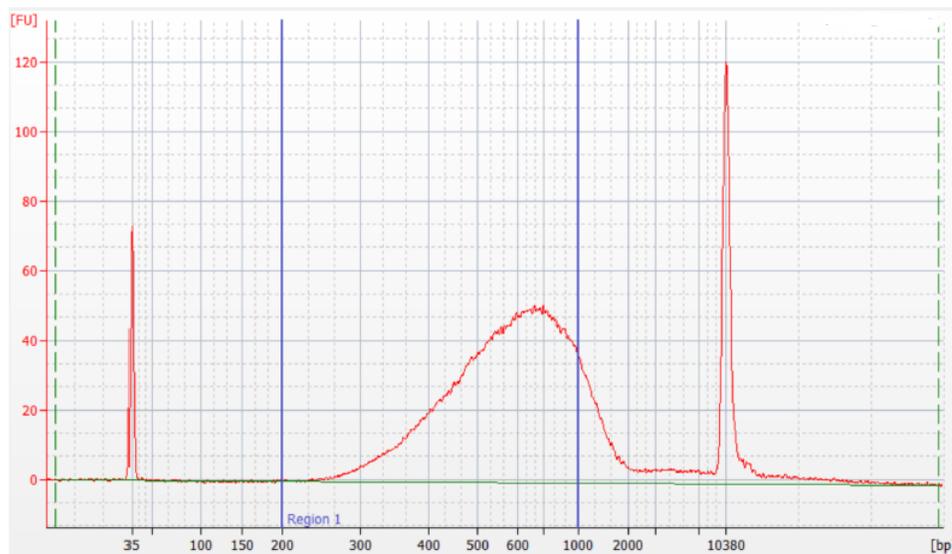

## 3. BD® AbSeq/Sample Tag library amplification

### 3.1 BD® AbSeq/Sample Tag PCR1

This section describes how to amplify BD® AbSeq/Sample Tag products through PCR.

1. In the pre-amplification workspace, pipet the following reagents into a new 1.5-mL LoBind® tube on ice.

#### BD® AbSeq/Sample Tag PCR1 Mix

Cap	Component	For 1 library (μL)	For 1 library with 20% overage (μL)	For 4 libraries with 20% overage (μL)	For 8 libraries with 20% overage (μL)
○	PCR MasterMix	100.0	120.0	480.0	960.0
○	Universal Oligo	10.0	12.0	48.0	96.0
●	Sample Tag PCR1 Primer	1.0	1.2	4.8	9.6
●	BD® AbSeq Primer	10.0	12.0	48.0	96.0
○	Nuclease-free water	12.0	14.4	57.6	115.2
	<b>Total</b>	<b>133.0</b>	<b>159.6</b>	<b>638.4</b>	<b>1276.8</b>

2. Gently vortex mix, briefly centrifuge, and place back on ice.

- In a new 1.5-mL tube, pipet 133  $\mu$ L of the BD<sup>®</sup> AbSeq/Sample Tag PCR1 Mix. Add 67  $\mu$ L of the BD<sup>®</sup> AbSeq/Sample Tag product from step 5 in [1.2 Denaturation and self-hybridization](#). Pipet-mix 10 times to create the BD<sup>®</sup> AbSeq/Sample Tag PCR1 Reaction Mix. Do not vortex.
- Pipet 50  $\mu$ L BD<sup>®</sup> AbSeq/Sample Tag PCR1 Reaction Mix into each of four 0.2-mL PCR tubes. Transfer any residual mix to one of the tubes.
- Bring the BD<sup>®</sup> AbSeq/Sample Tag PCR1 Reaction Mix to the post-amplification workspace.
- Run the following PCR program on the thermal cycler.

#### PCR1 conditions for BD<sup>®</sup> AbSeq/Sample Tag panel

Step	Cycles	Temperature	Time
Hot start	1	95 °C	3 min
Denaturation	10-11 <sup>a</sup>	95 °C	30 seconds
Annealing		60 °C	30 seconds
Extension		72 °C	1 min
Final extension	1	72 °C	5 min
Hold	1	4 °C	$\infty$

a. Recommended PCR cycles might need to be optimized for different cell types and cell number.

#### Recommended number of PCR cycles

Number of cells in PCR1	Recommended PCR cycles
7,500–10,000	11
20,000	10

 **STOPPING POINT:** The PCR can run overnight.

- After PCR, briefly centrifuge the tubes.
- Pipet-mix and combine the four reactions into a new 1.5-mL LoBind<sup>®</sup> tube, labeled AbSeq/Sample Tag PCR1. Keep the tube on ice.

### 3.2 BD<sup>®</sup> AbSeq/Sample Tag PCR1 cleanup

This section describes how to perform a single-sided AMPure<sup>®</sup> cleanup to remove primer dimers from the BD<sup>®</sup> AbSeq/Sample Tag PCR1 products. The final product is purified double-stranded DNA.

**Note:** Perform the purification in the post-amplification workspace.

- In a new 5.0-mL LoBind<sup>®</sup> tube, prepare 5 mL of fresh 80% (v/v) ethyl alcohol by combining 4 mL absolute ethyl alcohol, molecular biology grade, with 1 mL nuclease-free water. Vortex the tube for 10 seconds to mix.

**Note:** Make fresh 80% ethyl alcohol and use it within 24 hours. The 80% ethyl alcohol volume should be adjusted depending on the number of libraries.


2. Bring the AMPure® XP magnetic beads to room temperature. Vortex on high speed for 1 minute until the beads are fully resuspended.
3. To 200 µL of BD® AbSeq/Sample Tag PCR1 products, pipet 280 µL AMPure® beads (from step 8 in [3.1 BD® AbSeq/Sample Tag PCR1](#)).
4. Pipet-mix 10 times. Incubate at room temperature for 5 minutes.
5. Place the 1.5-mL LoBind® tube on the magnet for 5 minutes. Discard the supernatant.
6. Keeping the tube on the magnet, gently add 500 µL of fresh 80% ethyl alcohol into the tube and incubate for 30 seconds. Discard the supernatant.
7. Repeat step 6 once for a total of two washes.
8. Keeping the tube on the magnet, use a small-volume pipette to remove and discard any residual supernatant from the tube.
9. Air-dry the beads at room temperature for 5 minutes.
10. Remove the tube from the magnet and resuspend the bead pellet in 50 µL of Elution Buffer. Vigorously pipet-mix until the beads are uniformly dispersed. Small clumps do not affect performance.
11. Incubate at room temperature for 2 minutes and briefly centrifuge.
12. Place the tube on the magnet until the solution is clear, usually ~30 seconds.
13. Pipet the eluate (~50 µL) into a new 1.5-mL LoBind® tube (purified BD® AbSeq/Sample Tag PCR1 products).

 **STOPPING POINT:** Store at 2–8 °C before proceeding within 24 hours or at –25 °C to –15 °C for up to 6 months.

### 3.3 BD® AbSeq/Sample Tag PCR1 quality check

1. Measure the yield of the largest peak of the BD® AbSeq/Sample Tag PCR1 product (~150 bp) by using the Agilent 2100 Bioanalyzer with the High Sensitivity Kit assay. Follow the manufacturer's instructions.
2. Based on the yield of the largest peak, dilute an aliquot of BD® AbSeq PCR1 product to 0.1–1.1 ng/µL with nuclease-free water before index PCR of BD® AbSeq PCR1 products. See [3.6 BD® AbSeq/Sample Tag index PCR](#).

**Note:** Sample Tag PCR1 product does not need to be diluted before Sample Tag PCR2.

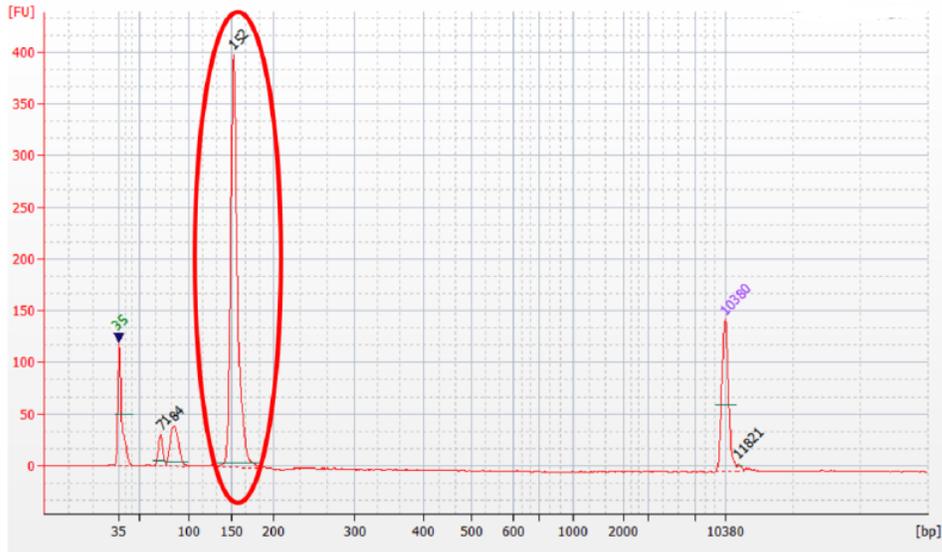**Figure 3** Sample Bioanalyzer High Sensitivity DNA trace - BD® AbSeq/Sample Tag PCR1

	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
▶ ◀	35	125.00	5,411.3	Lower Marker
2	71	51.24	1,087.6	
3	84	123.99	2,225.8	
4	152	838.65	8,334.2	
5	10,380	75.00	10.9	Upper Marker
6	11,821	0.00	0.0	

### 3.4 Sample Tag PCR2

This section describes how to amplify Sample Tag products through PCR. The PCR primers include partial Illumina sequencing adapters that enable the additions of full-length Illumina sequencing indices in the next PCR.

1. In the pre-amplification workspace, pipet reagents into a new 1.5-mL LoBind® tube on ice.

#### Sample Tag PCR2 Mix

Cap	Component	For 1 library (μL)	For 1 library with 20% overage (μL)	For 4 libraries with 20% overage (μL)	For 8 libraries with 20% overage (μL)
○	PCR MasterMix	25.0	30.0	120.0	240.0
○	Universal Oligo	2.0	2.4	9.6	19.2
●	Sample Tag PCR2 Primer	3.0	3.6	14.4	28.8
○	Nuclease-free water	15.0	18.0	72.0	144.0
	<b>Total</b>	<b>45.0</b>	<b>54.0</b>	<b>216.0</b>	<b>432.0</b>

2. Gently vortex mix, briefly centrifuge, and place back on ice. Pipet 45 μL of Sample Tag PCR2 Mix into one 0.2-mL PCR tube for each library.
3. Bring the Sample Tag PCR2 Mix to the post-amplification workspace.
4. Pipet 5 μL of PCR1 products into 45 μL Sample Tag PCR2 Mix to create the Sample Tag PCR2 Reaction Mix.
5. Gently vortex and briefly centrifuge.
6. Run the following PCR program on the thermal cycler.

#### PCR2 conditions for Sample Tag

Step	Cycles	Temperature	Time
Hot start	1	95 °C	3 min
Denaturation	10 <sup>a</sup>	95 °C	30 seconds
Annealing		66 °C	30 seconds
Extension		72 °C	1 min
Final extension	1	72 °C	5 min
Hold	1	4 °C	∞

a. Recommended PCR cycles might need to be optimized for different cell types and cell number.

 **STOPPING POINT:** The PCR can run overnight.

### 3.5 Sample Tag PCR2 cleanup and quantify

This section describes how to perform a single-sided AMPure® cleanup to remove primer dimers from the Sample Tag PCR2 products. The final product is purified double-stranded DNA.

**Note:** Perform PCR2 purification in the post-amplification workspace.

1. In a new 5.0-mL LoBind<sup>®</sup> tube, prepare 5 mL fresh 80% (v/v) ethyl alcohol by combining 4 mL absolute ethyl alcohol, molecular biology grade, with 1 mL of nuclease-free water. Vortex the tube for 10 seconds to mix.  
**Note:** Make fresh 80% ethyl alcohol and use it within 24 hours. The 80% ethyl alcohol volume should be adjusted depending on the number of libraries.
2. Bring AMPure<sup>®</sup> XP beads to room temperature and vortex at high speed for 1 minute until beads are fully resuspended.
3. Pipet 60  $\mu$ L of AMPure<sup>®</sup> beads into the Sample Tag PCR2 product.
4. Pipet-mix 10 times and incubate at room temperature for 5 minutes.
5. Place the tube on the strip tube magnet for 3 minutes. Discard the supernatant.
6. Keeping the tube on the magnet, gently add 200  $\mu$ L of fresh 80% ethyl alcohol into the tube and incubate for 30 seconds. Discard the supernatant.
7. Repeat step 6 once for a total of two washes.
8. Keeping the tube on the magnet, use a small-volume pipette to remove and discard any residual supernatant from the tube.
9. Air-dry the beads at room temperature for 3 minutes.
10. Remove the tube from the magnet and resuspend the bead pellet in 50  $\mu$ L of Elution Buffer. Pipet-mix until the beads are fully resuspended.
11. Incubate at room temperature for 2 minutes and briefly centrifuge.
12. Place the tube on the magnet until the solution is clear, usually ~30 seconds.
13. Pipet the eluate (~50  $\mu$ L) into a new 1.5-mL LoBind<sup>®</sup> tube.

 **STOPPING POINT:** Store at 2–8 °C before proceeding on the same day, or at –25 °C to –15 °C for up to 6 months.

14. Estimate the concentration by quantifying 2  $\mu$ L of the final sequencing library with a Qubit<sup>™</sup> Fluorometer using the Qubit<sup>™</sup> dsDNA HS Assay Kit. Follow the manufacturer's instructions.

### 3.6 BD<sup>®</sup> AbSeq/Sample Tag index PCR

This section describes how to generate BD<sup>®</sup> AbSeq/Sample Tag libraries compatible with the Illumina sequencing platform, by adding full-length Illumina sequencing adapters and indices through PCR.

1. In the pre-amplification workspace, pipet reagents into a new 1.5-mL LoBind® tube on ice.

#### BD® AbSeq/Sample Tag Index PCR Mix

Cap	Component	For 1 library (µL)	For 1 library with 20% overage (µL)	For 4 libraries with 20% overage (µL)	For 8 libraries with 20% overage (µL)
○	PCR MasterMix	25.0	30.0	120.0	240.0
●	Library Forward Primer	2.0	2.4	9.6	19.2
●	Library Reverse Primer 1–4 <sup>a</sup>	2.0	2.4	–	–
○	Nuclease-free water	18.0	21.6	86.4	172.8
	<b>Total</b>	<b>47.0</b>	<b>56.4</b>	<b>216.0</b>	<b>432.0</b>

a. For more than one library, use different Library Reverse Primers for each library.

2. Gently vortex mix, briefly centrifuge, and place back on ice.
3. Bring the BD® AbSeq/Sample Tag Index PCR Mix to post-amplification workspace.
4. In new 0.2 mL PCR tubes,
  - For BD® AbSeq library, pipet 3.0 µL of 0.1–1.1 ng/µL BD® AbSeq PCR1 product into 47 µL BD® AbSeq Index PCR Mix (from [3.3 BD® AbSeq/Sample Tag PCR1 quality check](#)). See the Concentration Index PCR input for BD® AbSeq/Sample Tag libraries (ng/µL) in [Recommended number of PCR cycles](#).
  - For Sample Tag libraries, pipet 3.0 µL of 0.1–1.1 ng/µL Sample Tag PCR2 products into 47.0 µL Sample Tag Index PCR Mix (from [3.5 Sample Tag PCR2 cleanup and quantify](#)). See the Concentration Index PCR input for BD® AbSeq/Sample Tag libraries (ng/µL) in [Recommended number of PCR cycles](#).
5. Gently vortex, and briefly centrifuge.
6. Run the following PCR program on the thermal cycler.

#### Index PCR conditions for BD® AbSeq/Sample Tag

Step	Cycles	Temperature	Time
Hot start	1	95 °C	3 min
Denaturation	Refer to the Recommended number of PCR cycles table <sup>a</sup>	95 °C	30 seconds
Annealing		60 °C	30 seconds
Extension		72 °C	30 seconds
Final extension	1	72 °C	1 min
Hold	1	4 °C	∞

a. Cycle number varies based on the concentration of the Sample Tag PCR2 product.

**Recommended number of PCR cycles**

Concentration Index PCR input for BD® AbSeq/Sample Tag libraries (ng/μL)	Recommended number of PCR cycles
0.5–1.1	6
0.25–0.5	7
0.1–0.25	8

 **STOPPING POINT:** The PCR can run overnight.

**3.7 BD® AbSeq/Sample Tag index PCR cleanup and quality check**


This section describes how to perform a single-sided AMPure® cleanup to remove primer dimers from the BD® AbSeq/Sample Tag Index PCR products. The final product is purified double-stranded DNA with full-length Illumina adapter sequences.

**Note:** Perform Index PCR purification in the post-amplification workspace.


- In a new 5.0-mL LoBind® tube, prepare 5 mL fresh 80% (v/v) ethyl alcohol by combining 4 mL absolute ethyl alcohol, molecular biology grade, with 1 mL of nuclease-free water. Vortex the tube for 10 seconds to mix.  
**Note:** Make fresh 80% ethyl alcohol, and use it within 24 hours. The 80% ethyl alcohol volume should be adjusted depending on the number of libraries.
- Bring AMPure® XP beads to room temperature and vortex at high speed for 1 minute until the beads are fully resuspended.
- Briefly centrifuge all the index PCR products.
- To 50.0 μL of the BD® AbSeq/Sample Tag index PCR products, pipet 40 μL AMPure® beads.
- Pipet-mix 10 times and incubate at room temperature for 5 minutes.
- Place the tube on the strip tube magnet for 3 minutes. Discard the supernatant.
- Keeping the tube on the magnet, gently add 200 μL of fresh 80% ethyl alcohol into the tube, and incubate for 30 seconds. Discard the supernatant.
- Repeat step 7 for a total of two washes.
- Keeping the tube on the magnet, use a small-volume pipette to remove and discard the residual supernatant from the tube.
- Air-dry the beads at room temperature for 3 minutes.
- Remove the tube from the magnet and resuspend the bead pellet in 50 μL of Elution Buffer. Pipet-mix until the beads are fully resuspended.
- Incubate at room temperature for 2 minutes, and briefly centrifuge.
- Place the tube on the magnet until the solution is clear, usually ~30 seconds.
- Pipet the entire eluate (~50 μL) into a new 1.5-mL LoBind® tube (final sequencing libraries).
- Perform quality control before freezing samples.

- a. Estimate the concentration by quantifying 2  $\mu\text{L}$  of the final sequencing library with a Qubit™ Fluorometer using the Qubit™ dsDNA HS Kit to obtain an approximate concentration of PCR products to dilute for quantification on an Agilent 2100 Bioanalyzer. Follow the manufacturer's instructions. The expected concentration of the libraries is  $>1.5 \text{ ng}/\mu\text{L}$ .
- b. Measure the average fragment size of the BD® AbSeq/Sample Tag library within the size range of 200–1,000 bp by using the Agilent Bioanalyzer with the High Sensitivity Kit for 50–7,000 bp, 5–1,000  $\text{pg}/\mu\text{L}$ . Follow the manufacturer's instructions.

**Figure 4** Sample Bioanalyzer High Sensitivity DNA trace - Sample Tag index PCR product

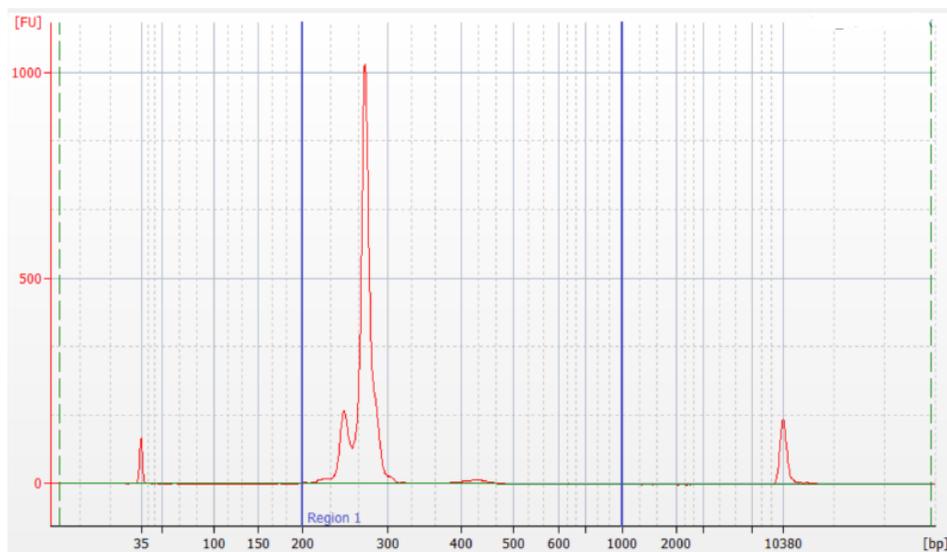

**Figure 5** Sample Bioanalyzer High Sensitivity DNA trace - BD® AbSeq index PCR product

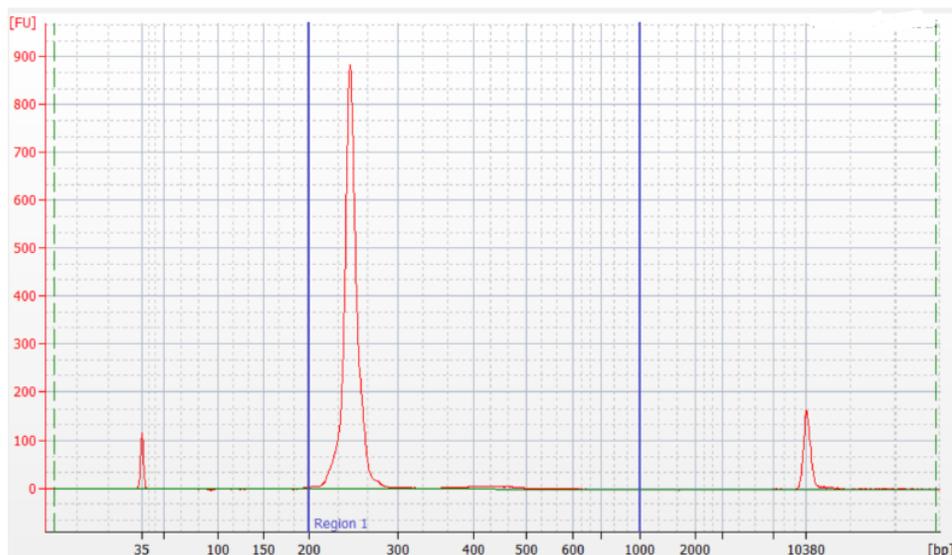

**STOPPING POINT:** Store at  $-25\text{ }^{\circ}\text{C}$  to  $-15\text{ }^{\circ}\text{C}$  for up to 6 months until sequencing.

## 4. TCR/BCR library amplification

### Before you begin

- Obtain Exonuclease I-treated BD Rhapsody™ Enhanced Cell Capture Beads.
- Thaw reagents in the BD Rhapsody™ Mouse TCR/BCR Next Amplification Kit at room temperature, and then place on ice.

### 4.1 TCR/BCR PCR1

1. Obtain beads from Step 19 of [2.1 Random priming and extension \(RPE\)](#).

**Note:** Use the entire sample of beads. Sub-sampling beads is not recommended for TCR/BCR FL + WTA + BD® AbSeq + Sample Tag combination assays.

2. In the pre-amplification workspace, pipet the following reagents into a new 1.5-mL LoBind® tube.

#### TCR/BCR PCR1 Mix

Cap	Component	For 1 library (μL)	For 1 library with 20% overage (μL)	For 4 libraries with 20% overage (μL)	For 8 libraries with 20% overage (μL)
○	PCR MasterMix	100.00	120.00	480.00	960.00
●	TCR/BCR Universal Oligo N1	10.00	12.00	48.00	96.00
●	Bead RT/PCR Enhancer	12.00	14.40	57.6	115.20
●	Mouse TCR N1 Primer <sup>a</sup>	2.40	2.88	11.52	23.04
●	Mouse BCR N1 Primer <sup>a</sup>	2.40	2.88	11.52	23.04
○	Nuclease-free water	73.20	87.84	351.36	702.72
	<b>Total</b>	<b>200.00</b>	<b>240.00</b>	<b>960.00</b>	<b>1920.00</b>

a. If only doing TCR or BCR amplification, replace N1 primer volume with nuclease-free water. For example, if only doing TCR amplification, replace BCR N1 primer with nuclease-free water.

3. Gently vortex mix, briefly centrifuge, and place back on ice.
4. Briefly spin the tube with the bead suspension. Place the tube of beads in Bead Resuspension Buffer on a 1.5-mL magnet for ≤1 minute. Discard the supernatant.
5. Remove the tube from the magnet and resuspend the beads in 200 μL of TCR/BCR PCR1 Mix to create the TCR/BCR PCR1 Reaction Mix. Do not vortex.
6. Ensuring that the beads are fully resuspended, pipet 50 μL of TCR/BCR PCR1 Reaction Mix with beads into each of four 0.2-mL PCR tubes. Transfer any residual mix to one of the tubes.
7. Bring the TCR/BCR PCR1 Reaction Mix to the post-amplification workspace.

8. Run the following PCR program on the thermal cycler.

**PCR1 conditions for TCR/BCR panel**

Step	Cycles	Temperature	Time
Hot start	1	95 °C <sup>a</sup>	3 min
Denaturation	10-11 <sup>b</sup>	95 °C	30 seconds
Annealing		60 °C	1 min
Extension		72 °C	1 min
Final extension	1	72 °C	5 min
Hold	1	4 °C	∞

a. To avoid beads settling due to prolonged incubation time on the thermal cycler before the denaturation step, it is critical to pause the instrument at 95 °C before loading the samples. Different thermal cyclers might have different pause time settings. In certain brands of thermal cyclers, however, we have observed a step-skipping error with the pause/unpause functions. To ensure that the full 3-minute denaturation is not skipped, verify that the pause/unpause functions are working correctly on your thermal cycler. To avoid the step-skipping problem, a 1-minute 95 °C pause step can be added immediately before the 3-minute 95 °C denaturation step.

b. Recommended PCR cycles might need to be optimized for different cell types and cell number.

**Recommended number of PCR cycles**

Number of cells in PCR1	Recommended PCR cycles
7,500 - 10,000	11
20,000	10

9. Ramp the heated lid and heat block of the post-amplification thermal cycler to  $\geq 95$  °C by starting the thermal cycler program and then pausing it.

**Note:** Do not proceed to thermal cycling until each tube is gently mixed by pipette to ensure uniform bead suspension.

10. For each 0.2-mL PCR tube, gently pipet-mix, immediately place the tube in thermal cycler, and unpause the thermal cycler program.

 **STOPPING POINT:** The PCR can run overnight, but proceed with purification within 24 hours after PCR.

11. After PCR, briefly centrifuge the tubes.
12. Put the tubes on a strip tube magnet for >30 seconds. Remove and carefully combine the supernatant in a new 1.5-mL LoBind<sup>®</sup> tube without disturbing the beads.

 **Important:** Save the supernatant for purification of TCR/BCR PCR1 products. Keep on ice until proceeding to TCR/BCR PCR2.

13. Place the 0.2-mL PCR tubes on a magnet for  $\leq 1$  minute. Retain the supernatant.
14. For each sample, collect and combine the supernatant from the four 0.2-mL PCR tubes into one new 1.5-mL LoBind<sup>®</sup> tube. Discard the beads.

## 4.2 TCR/BCR PCR1 cleanup

This section describes how to perform a single-sided AMPure<sup>®</sup> cleanup to remove primer dimers from the TCR/BCR PCR1 products. The final product is purified double-stranded DNA.

**Note:** Perform the purification in the post-amplification workspace.

1. In a new 5.0-mL LoBind<sup>®</sup> tube, prepare 5 mL of fresh 80% (v/v) ethyl alcohol by combining 4 mL absolute ethyl alcohol, molecular biology grade, with 1 mL nuclease-free water. Vortex the tube for 10 seconds to mix.  
**Note:** Make fresh 80% ethyl alcohol and use it within 24 hours. The 80% ethyl alcohol volume should be adjusted depending on the number of libraries.
2. Bring the AMPure<sup>®</sup> XP magnetic beads to room temperature. Vortex on high speed for 1 minute until the beads are fully resuspended.
3. To 200  $\mu$ L of TCR/BCR PCR1 products, pipet 140  $\mu$ L AMPure<sup>®</sup> beads (from step 13 in [4.1 TCR/BCR PCR1](#)).
4. Pipet-mix 10 times. Incubate at room temperature for 5 minutes.
5. Place the 1.5-mL LoBind<sup>®</sup> tube on the magnet for 5 minutes. Discard the supernatant.
6. Keeping the tube on the magnet, gently add 500  $\mu$ L of fresh 80% ethyl alcohol into the tube and incubate for 30 seconds. Discard the supernatant.
7. Repeat step 6 once for a total of two washes.
8. Keeping the tube on the magnet, use a small-volume pipette to remove and discard any residual supernatant from the tube.
9. Air-dry the beads at room temperature for 5 minutes.
10. Remove the tube from the magnet and resuspend the bead pellet in 50  $\mu$ L of Elution Buffer. Vigorously pipet-mix until the beads are uniformly dispersed. Small clumps do not affect performance.
11. Incubate at room temperature for 2 minutes and briefly centrifuge.
12. Place the tube on the magnet until the solution is clear, usually ~30 seconds.
13. Pipet the eluate (~50  $\mu$ L) into a new 1.5-mL LoBind<sup>®</sup> tube (purified TCR/BCR PCR1 products).

 **STOPPING POINT:** Store at 2–8 °C before proceeding within 24 hours or at –25 °C to –15 °C for up to 6 months.

## 4.3 TCR/BCR PCR2

This section describes how to amplify TCR/BCR products through PCR.

1. In the pre-amplification workspace, pipet reagents into a new 1.5-mL LoBind® tube on ice.

#### TCR and BCR PCR2 Mixes

Cap	Component	For 1 library (μL)	For 1 library with 20% overage (μL)	For 4 libraries with 20% overage (μL)	For 8 libraries with 20% overage (μL)
○	PCR MasterMix	25.0	30.0	120.0	240.0
●	TCR/BCR Universal Oligo N2	2.0	2.4	9.6	19.2
● ●	Mouse TCR N2 Primer or Mouse BCR N2 Primer <sup>a</sup>	6.0	7.2	28.8	57.6
○	Nuclease-free water	12.0	14.4	57.6	115.2
	<b>Total</b>	<b>45.0</b>	<b>54.0</b>	<b>216.0</b>	<b>432.0</b>

a. PCR2 mixes for TCR and BCR are made separately.

2. Gently vortex mix, briefly centrifuge, and place back on ice. Pipet 45 μL of PCR2 Mix into one 0.2-mL PCR tube for each library.
3. Bring the TCR PCR2 Mix and the BCR PCR2 Mix to the post-amplification workspace.
4. Pipet 5.0 μL of PCR1 products into 45 μL of PCR2 Mix for each library to create the TCR PCR2 Reaction Mix and BCR PCR2 Reaction Mix, respectively.
5. Gently vortex and briefly centrifuge.
6. Run the following PCR program on the thermal cycler.

**PCR2 conditions for TCR and BCR**

Step	Cycles	Temperature	Time	
Phase I:	1	95 °C	3 min	Press Option > <b>Auto Delta Starting cycle</b> > "2" > Delta > "1 degree" > Done The temperature decreases by 1 °C each cycle, from 70 °C to 56 °C.
	15	95 °C	30 seconds	
		70–56 °C	1 min	
		72 °C	1 min	
Phase II:	8	95 °C	30 seconds	
		55 °C	1 min	
		72 °C	1 min	
	1	72 °C	5 min	
	1	4 °C	∞	

 **STOPPING POINT:** The PCR can run overnight.

**4.4 TCR/BCR PCR2 cleanup**

This section describes how to perform a single-sided AMPure® cleanup to remove primer dimers from the TCR and BCR PCR2 products. The final product is purified double-stranded DNA.

**Note:** Perform PCR2 purification in the post-amplification workspace.

1. In a new 5.0-mL LoBind® tube, prepare 5 mL fresh 80% (v/v) ethyl alcohol by combining 4 mL absolute ethyl alcohol, molecular biology grade, with 1 mL of nuclease-free water. Vortex the tube for 10 seconds to mix.

**Note:** Make fresh 80% ethyl alcohol and use it within 24 hours. The 80% ethyl alcohol volume should be adjusted depending on the number of libraries.

2. Bring AMPure® XP beads to room temperature and vortex at high speed for 1 minute until beads are fully resuspended.
3. To 50.0 µL PCR2 products, pipet 35 µL of AMPure® beads.
4. Pipet-mix 10 times and incubate at room temperature for 5 minutes.
5. Place the tube on the strip tube magnet for 3 minutes. Discard the supernatant.
6. Keeping the tube on the magnet, gently add 200 µL of fresh 80% ethyl alcohol into the tube and incubate for 30 seconds. Discard the supernatant.
7. Repeat step 6 once for a total of two washes.
8. Keeping the tube on the magnet, use a small-volume pipette to remove and discard any residual supernatant from the tube.

9. Air-dry the beads at room temperature for 3 minutes.
10. Remove the tube from the magnet and resuspend the bead pellet in 50  $\mu\text{L}$  of Elution Buffer. Pipet-mix until the beads are fully resuspended.
11. Incubate at room temperature for 2 minutes and briefly centrifuge.
12. Place the tube on the magnet until the solution is clear, usually ~30 seconds.
13. Pipet the eluate (~50  $\mu\text{L}$ ) into a new 1.5-mL LoBind® tube.

**⚠ STOPPING POINT:** Store at 2–8 °C before proceeding on the same day, or at –25 °C to –15 °C for up to 6 months.

14. Estimate the concentration by quantifying 2  $\mu\text{L}$  of the final sequencing library with a Qubit™ Fluorometer using the Qubit™ dsDNA HS Assay Kit. Follow the manufacturer's instructions.

#### 4.5 TCR/BCR RPE

1. Dilute an aliquot of the TCR and BCR PCR2 products with nuclease-free water to 1.0 ng/ $\mu\text{L}$ .

**Note:** If PCR2 concentration is <1 ng/ $\mu\text{L}$ , increase the volume of PCR2 product needed to ensure 5 ng total concentration and decrease the volume of water in the Random Primer Mix accordingly.

2. In pre-amplification workspace, pipet the following reagents into a new 1.5 mL LoBind® tube:

##### Random Primer Mix

Cap	Component	For 1 library ( $\mu\text{L}$ )	For 1 library with 20% overage ( $\mu\text{L}$ )	For 4 libraries with 20% overage ( $\mu\text{L}$ )	For 8 libraries with 20% overage ( $\mu\text{L}$ )
○	TCR/BCR Extension Buffer	5.0	6.0	24.0	48.0
○	TCR/BCR Extension Primers	2.5	3.0	12.0	24.0
○	Nuclease-free water	Up to 34.0	Up to 40.8	Up to 163.2	Up to 326.4
	<b>Total</b>	<b>41.5</b>	<b>49.8</b>	<b>199.2</b>	<b>398.4</b>

3. Pipet-mix the Random Primer Mix and keep at room temperature. Pipet 41.5  $\mu\text{L}$  of Random Primer Mix into one 0.2-mL PCR tube for each library.
4. Add 5  $\mu\text{L}$  of 1.0 ng/ $\mu\text{L}$  purified TCR or BCR PCR2 products into each 0.2-mL PCR tube containing Random Primer Mix.
5. Perform denaturation and random priming on thermocycler using the following program:

**Program**

Temperature	Time	Cycles
95 °C	5 min	1
37 °C	5 min	
25 °C	15 min	

- Briefly centrifuge the tube and keep at room temperature.
- In pre-amplification workspace, pipet the following reagents into a new 1.5 mL LoBind<sup>®</sup> tube:

**Primer Extension Enzyme Mix**

Cap	Component	For 1 library (μL)	For 1 library with 20% overage (μL)	For 4 libraries with 20% overage (μL)	For 8 libraries with 20% overage (μL)
	dNTP	2.0	2.4	9.6	19.2
	TCR/BCR Extension Enzyme	1.5	1.8	7.2	14.4
	<b>Total</b>	<b>3.5</b>	<b>4.2</b>	<b>16.8</b>	<b>33.6</b>

- Gently vortex mix, centrifuge, and place on ice.
- Add 3.5 μL Primer Extension Enzyme Mix to the Random Priming Reaction tube to bring total volume up to 50 μL. Run the following protocol on a thermocycler for Extension:

**Program**

Temperature	Time	Cycles
25 °C	10 min	1
37 °C	15 min	
45 °C	10 min	
55 °C	10 min	

- Remove the tubes from the thermocycler and prepare to purify RPE product.

**4.6 TCR/BCR RPE cleanup**

**Note:** Perform purification in the post-amplification workspace.

- In a new 5.0-mL LoBind<sup>®</sup> tube, prepare 5 mL fresh 80% (v/v) ethyl alcohol by combining 4 mL absolute ethyl alcohol, molecular biology grade, with 1 mL of nuclease-free water. Vortex the tube for 10 seconds to mix.

**Note:** Make fresh 80% ethyl alcohol and use it within 24 hours. The 80% ethyl alcohol volume should be adjusted depending on the number of libraries.

2. Bring AMPure® XP beads to room temperature and vortex at high speed for 1 minute until beads are fully resuspended.
3. To the TCR and BCR RPE products, add 90 µL AMPure® beads.
4. Pipet-mix 10 times and incubate at room temperature for 5 minutes.
5. Place the tube on the strip tube magnet for 3 minutes. Discard the supernatant.
6. Keeping the tube on the magnet, gently add 200 µL of fresh 80% ethyl alcohol into the tube and incubate for 30 seconds. Discard the supernatant.
7. Repeat step 6 once for a total of two washes.
8. Keeping the tube on the magnet, use a small-volume pipette to remove and discard any residual supernatant from the tube.
9. Air-dry the beads at room temperature for 1 minute.
10. Remove tubes from the magnet and add 50 µL of Elution Buffer.
11. Incubate at room temperature for 2 minutes and briefly centrifuge.
12. Place the tube on the magnet until the solution is clear, usually ≤30 seconds.
13. Pipet the entire eluate (~50 µL) into a new 1.5-mL LoBind® tube separately (purified TCR and BCR RPE products).

#### 4.7 TCR/BCR index PCR

This section describes how to generate TCR/BCR libraries compatible with the Illumina sequencing platform, by adding full-length Illumina sequencing adapters and indices through PCR.

1. In the pre-amplification workspace, pipet reagents into a new 1.5-mL LoBind® tube on ice.

##### TCR/BCR Index PCR Mix

Cap	Component	For 1 library (µL)	For 1 library with 20% overage (µL)	For 4 libraries with 20% overage (µL)	For 8 libraries with 20% overage (µL)
○	PCR MasterMix	25.0	30.0	120.0	240.0
●	Library Forward Primer <sup>a</sup>	2.0	2.4	9.6	19.2
●	Library Reverse Primer <sup>a</sup>	2.0	2.4	–	–
	<b>Total</b>	<b>29.0</b>	<b>34.8</b>	<b>129.6</b>	<b>259.2</b>

a. For more than one library, use different Library Reverse Primers for each TCR or BCR library. If needed, use the BD Rhapsody™ Dual Index Kit to process additional samples.

2. Gently vortex mix, briefly centrifuge, and place back on ice.
3. Bring the TCR/BCR Index PCR Mix to post-amplification workspace.

4. In new 0.2 mL PCR tubes, pipet 21  $\mu\text{L}$  of TCR/BCR RPE purified products into 29  $\mu\text{L}$  of TCR/BCR Index PCR Mix.
5. Gently vortex, and briefly centrifuge.
6. Run the following PCR program on the thermal cycler.

#### Index PCR conditions for TCR/BCR

Step	Cycles	Temperature	Time
Hot start	1	95 °C	3 min
Denaturation	10	95 °C	30 seconds
Annealing		60 °C	30 seconds
Extension		72 °C	30 seconds
Final extension	1	72 °C	1 min
Hold	1	4 °C	$\infty$

 **STOPPING POINT:** The PCR can run overnight.

### 4.8 TCR/BCR index PCR cleanup and quality check

This section describes how to perform a single-sided AMPure<sup>®</sup> cleanup to remove primer dimers from the TCR/BCR index PCR products. The final product is purified double-stranded DNA with full-length Illumina adapter sequences.

**Note:** Perform Index PCR purification in the post-amplification workspace.


1. In a new 5.0-mL LoBind<sup>®</sup> tube, prepare 5 mL fresh 80% (v/v) ethyl alcohol by combining 4 mL absolute ethyl alcohol, molecular biology grade, with 1 mL of nuclease-free water. Vortex the tube for 10 seconds to mix.


**Note:** Make fresh 80% ethyl alcohol, and use it within 24 hours. The 80% ethyl alcohol volume should be adjusted depending on the number of libraries.

2. Bring AMPure<sup>®</sup> XP beads to room temperature and vortex at high speed for 1 minute until the beads are fully resuspended.
3. Briefly centrifuge all the index PCR products.
4. Transfer 40  $\mu\text{L}$  of the TCR and/or BCR index PCR products to a new strip tube or tubes.
5. Pipet 26  $\mu\text{L}$  of AMPure<sup>®</sup> beads.
6. Pipet-mix 10 times and incubate at room temperature for 5 minutes.
7. Place the tube on the strip tube magnet for 3 minutes. Discard the supernatant.
8. Keeping the tube on the magnet, gently add 200  $\mu\text{L}$  of fresh 80% ethyl alcohol into the tube and incubate for 30 seconds. Discard the supernatant.
9. Repeat step 7 for a total of two washes.

10. Keeping the tube on the magnet, use a small-volume pipette to remove and discard the residual supernatant from the tube.
11. Air-dry the beads at room temperature for 3 minutes.
12. Remove the tube from the magnet and resuspend the bead pellet in 50  $\mu\text{L}$  of Elution Buffer. Pipet-mix until the beads are fully resuspended.
13. Incubate at room temperature for 2 minutes, and briefly centrifuge.
14. Place the tube on the magnet until the solution is clear, usually  $\sim 30$  seconds.
15. Pipet the entire eluate ( $\sim 50 \mu\text{L}$ ) into a new 1.5-mL LoBind® tube (final sequencing libraries).
16. Perform quality control before freezing samples.
  - a. Estimate the concentration by quantifying 2  $\mu\text{L}$  of the final sequencing library with a Qubit™ Fluorometer using the Qubit™ dsDNA HS Kit to obtain an approximate concentration of PCR products to dilute for quantification on an Agilent 2100 Bioanalyzer. Follow the manufacturer's instructions. The expected concentration of the libraries is  $>1.5 \text{ ng}/\mu\text{L}$ .
  - b. Measure the average fragment size of the TCR/BCR library within the size range of 200–1,000 bp by using the Agilent Bioanalyzer with the High Sensitivity Kit for 50–7,000 bp, 5–1,000 pg/ $\mu\text{L}$ . Follow the manufacturer's instructions.

**Figure 6** Sample Bioanalyzer High Sensitivity DNA trace - TCR index PCR product

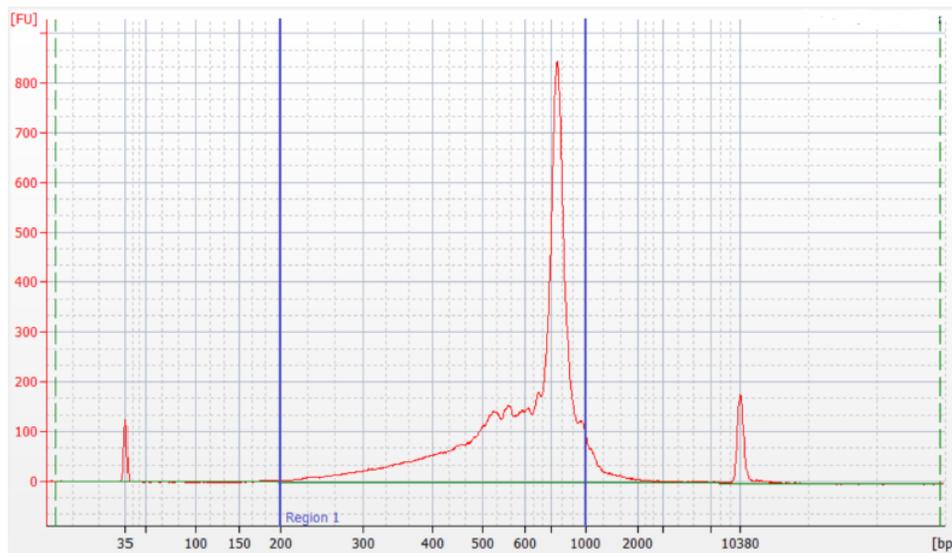

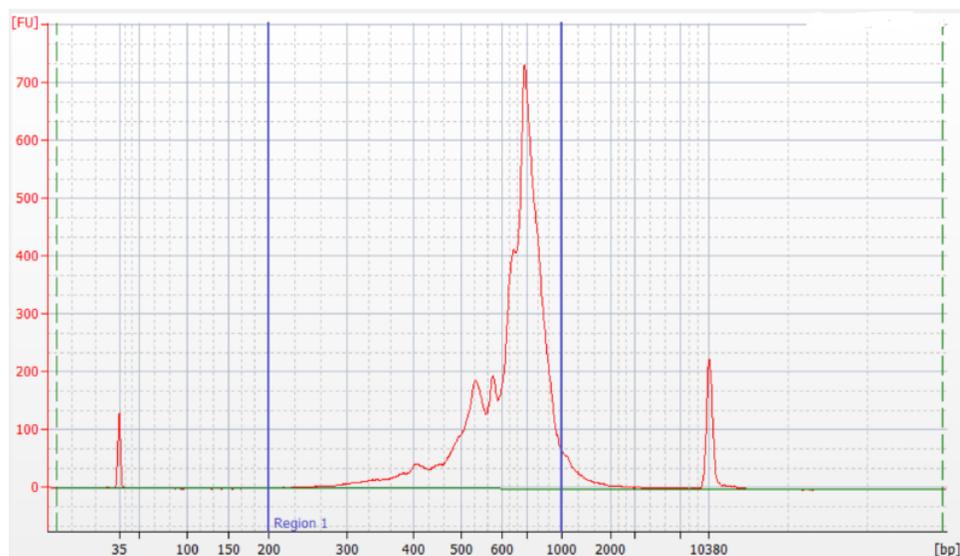**Figure 7** Sample Bioanalyzer High Sensitivity DNA trace - BCR index PCR product

**⚠ STOPPING POINT:** Store at  $-25\text{ }^{\circ}\text{C}$  to  $-15\text{ }^{\circ}\text{C}$  for up to 6 months until sequencing.

## Sequencing

The sequencing depth for each library is dependent on application. For cell-type clustering, shallow sequencing is sufficient. However, for in-depth analysis, such as comparison across multiple samples, deep sequencing is advised. We recommend meeting the requirement for recursive substitution error correction (RSEC) sequencing depth of  $\geq 6$  to reach the threshold of sequencing saturation where most molecules of the library have been recovered, approximately 80%. The RSEC sequencing depth and sequencing saturation are both reported by the analysis pipeline. The actual sequencing reads/cell required to achieve this depth can vary as it is dependent on the gene expression levels, number of cells, and sequencing run quality. The following reads/cell are recommended for WTA mRNA, Sample Tag, BD<sup>®</sup> AbSeq, TCR, and BCR libraries.

### Read requirements for libraries

Library	Read requirement for data analysis
WTA mRNA	10,000–100,000 reads/cell
Sample Tag: Samples from the same type of cell	120 reads/cell
Sample Tag: Samples from different types of cells	600 reads/cell
BD <sup>®</sup> AbSeq	1,000 reads/cell/AbSeq <sup>a</sup>
TCR	~5,000 reads/T cell
BCR	~5,000 reads/B cell

a. The amount of sequencing needed for BD<sup>®</sup> AbSeq libraries will vary depending on application, BD<sup>®</sup> AbSeq panel plexy, and cell type.

## Pooling libraries for sequencing

The efficiency of sequencing on Illumina instruments is influenced by many conditions, library size being one of them. The TCR and BCR libraries are ~200–300 bp larger than the WTA mRNA library and ~600 bp larger than the BD® AbSeq and Sample Tag libraries which will cause them to produce less sequencing data if pooled in a 1:1 ratio with the other libraries. To overcome the difference in sequencing efficiency, more DNA of the TCR and BCR libraries needs to be included in the pool than would be expected when calculating ratios based on read depth. The following tables show examples of different pooling strategies and the expected sequencing outcome, with and without correction for the size of the TCR and BCR libraries. Validation data indicates a 5X volume correction factor is needed for sequencing TCR and BCR libraries with WTA mRNA, BD® AbSeq, and Sample Tag libraries.

**⚠ Important:** BD® AbSeq and Sample Tag libraries can be sequenced together or separately from WTA and TCR/BCR libraries. For optimal clustering of TCR/BCR libraries on Illumina platforms, however, we recommend sequencing BD® AbSeq and Sample Tag libraries separately. Recommendations for sequencing BD® AbSeq and Sample Tag libraries alone or in a pool with other libraries are in the next section.

### Example of pooling with no correction

In this example, a total of 5,000 enriched T cells were processed as a combination of two samples (two sample tags) with a 10-plex BD® AbSeq panel. These calculations assume the TCR library, and BCR library if included, sequences at 1/5 the efficiency of the WTA mRNA, BD® AbSeq, and Sample Tag libraries, supported by internal testing.

#### Pooling for WTA mRNA, TCR, BD® AbSeq, and Sample Tag libraries

A	B	C	D	E	F	G
Library type	Number of cells	Expected reads/cell	Reads needed	Pooling ratio	Sequencing results	Sequencing results (reads/cell)
WTA mRNA	5,000	25,000	125,000,000	62.3%	139,236,902	27,847
TCR	5,000	5,000	25,000,000	12.5%	5,000,000	1,000
BD® AbSeq	5,000	10,000	50,000,000	24.9%	55,694,761	11,139
Sample Tag	5,000	120	600,000	0.3%	668,337	134
<b>Total</b>			<b>200,600,000</b>	<b>100%</b>	<b>200,600,000</b>	–

After sequencing, the reads/cell for the TCR library (Column G) does not match with and are much lower than the expectation (Column C), because it does not sequence as efficiently as the WTA mRNA, BD® AbSeq, nor Sample Tag library. The remaining reads are allotted to the WTA mRNA, BD® AbSeq, and Sample Tag libraries resulting in more reads than required. To obtain the desired number of reads/cell for each library, a correction factor is required for pooling calculations.

### Example of pooling with correction

In this example, the same sample as in the previous one was pooled using a correction factor of 5 for the TCR library to overcome the differences in sequencing efficiency. The amount of data needing to be generated (Column D) is based on the cell number (Column B) and expected number of reads per cell (Column C). Based on this example, 200.6 million reads are needed to achieve the appropriate read depths. Changing the pooling ratios by correcting for the lower TCR sequencing efficiency will help ensure the correct amount of data is

generated for each library. This modified pooling scheme, however, does not change the total amount of data needing to be generated, which is 200.6 million reads.

### Pooling for WTA mRNA, TCR, BD<sup>®</sup> AbSeq, and Sample Tag libraries

A	B	C	D	E	F	G	H	I	J
Library type	Number of cells	Expected reads/cell	Reads needed	Pooling ratio before correction	Correction	Reads needed for pooling	Pooling ratio with correction	Sequencing results	Sequencing results (reads/cell)
WTA mRNA	5,000	25,000	125,000,000	62.3%	N/A	125,000,000	41.6%	125,000,000	25,000
TCR	5,000	5,000	25,000,000	12.5%	5 <sup>a</sup>	125,000,000	41.6%	25,000,000	5,000
BD <sup>®</sup> AbSeq	5,000	10,000	50,000,000	24.9%	N/A	50,000,000	16.6%	50,000,000	10,000
Sample Tag	5,000	120	600,000	0.3%	N/A	600,000	0.2%	600,000	120
<b>Total</b>			<b>200,600,000<sup>b</sup></b>	<b>100%</b>	–	<b>300,600,000<sup>c</sup></b>	<b>100%</b>	<b>200,600,000</b>	–

a. The 5X correction factor is a recommended starting point and some fine tuning may be required to achieve the optimal library balance.

b. Total amount of data to be requested from the sequencing facility plus 3% PhiX.

c. Read total only for pooling purposes.

After sequencing, the total amount of data generated (Column I) as well as the reads/cell for each library (Column J) are as expected (Columns D and C). The correction for library pooling did not change the amount of data generated (200.6 million reads) but helped ensure the data was spread out appropriately to each library.

### Example of pooling with a mixed population

The following table shows the pooling logic for two samples of purified B and T cells, each stained with a different sample tag and 10-plex BD<sup>®</sup> AbSeq panel.

**Pooling for WTA mRNA, TCR, BCR, BD® AbSeq, and Sample Tag libraries**

A	B	C	D	E	F	G	H	I	J
Library type	Number of cells	Expected reads/cell	Reads needed	Pooling ratio before correction	Correction	Reads needed for pooling	Pooling ratio with correction	Sequencing results	Sequencing results (reads/cell)
WTA mRNA	10,000	25,000	250,000,000	62%	N/A	250,000,000	41%	250,000,000	25,000
TCR	7,000	5,000	35,000,000	9%	5 <sup>a</sup>	175,000,000	29%	35,000,000	5,000
BCR	3,000	5,000	15,000,000	3%	5 <sup>a</sup>	75,000,000	12%	15,000,000	5,000
BD® AbSeq	10,000	10,000	100,000,000	25%	N/A	100,000,000	17%	100,000,000	10,000
Sample Tag	10,000	600	6,000,000	1%	N/A	6,000,000	1%	6,000,000	600
<b>Total</b>			<b>406,000,000<sup>b</sup></b>	<b>100%</b>	–	<b>606,000,000<sup>c</sup></b>	<b>100%</b>	<b>406,000,000</b>	–

a. The 5X correction factor is a recommended starting point and some fine tuning may be required to achieve the optimal library balance.

b. Total amount of data to be requested from the sequencing facility plus 3% PhiX.

c. Read total only for pooling purposes.

**Example of pooling with correction while sequencing BD® AbSeq and Sample Tag libraries separately**

The following table shows the pooling logic for the previous example of a mixed population where the BD® AbSeq and Sample Tag libraries are sequenced together and the WTA mRNA, TCR, and BCR libraries are sequenced together.

**Pooling for BD® AbSeq and Sample Tag libraries**

A	B	C	D	E	F	G
Library type	Number of cells	Expected reads/cell	Reads needed	Pooling ratio	Sequencing results	Sequencing results (reads/cell)
BD® AbSeq	10,000	10,000	100,000,000	94%	100,000,000	10,000
Sample Tag	10,000	600	6,000,000	6%	6,000,000	600
<b>Total</b>			<b>106,000,000<sup>a</sup></b>	<b>100%</b>	<b>106,000,000</b>	–

a. Total amount of data to be requested from the sequencing facility plus 1% PhiX.

### Pooling for WTA mRNA, TCR, and BCR libraries

A	B	C	D	E	F	G	H	I	J
Library type	Number of cells	Expected reads/cell	Reads needed	Pooling ratio before correction	Correction	Reads needed for pooling	Pooling ratio with correction	Sequencing results	Sequencing results (reads/cell)
WTA mRNA	10,000	25,000	250,000,000	83%	N/A	250,000,000	50%	250,000,000	25,000
TCR	7,000	5,000	35,000,000	12%	5 <sup>a</sup>	175,000,000	35%	35,000,000	5,000
BCR	3,000	5,000	15,000,000	5%	5 <sup>a</sup>	75,000,000	15%	15,000,000	5,000
<b>Total</b>			<b>300,000,000<sup>b</sup></b>	<b>100%</b>	–	<b>500,000,000<sup>c</sup></b>	<b>100%</b>	<b>300,000,000</b>	–

a. The 5X correction factor is a recommended starting point and some fine tuning may be required to achieve the optimal library balance.

b. Total amount of data to be requested from the sequencing facility plus 3% PhiX.

c. Read total only for pooling purposes.

### Additional considerations

- The 5X volume correction factor is a recommended starting place for pooling these libraries. This may need to be adjusted to accommodate different types of flow cells (for example, patterned vs non-patterned).
- It can be easier to achieve the desired sequencing depth when sequencing multiple TCR or BCR libraries alone since all the libraries are the same size. Pooling will not require a correction and will only be dependent on the number of cells and the reads/cell. This scheme, however, would require 10–15% PhiX, rather than the 3% PhiX when sequencing with the WTA mRNA library.
- Sequencing BD<sup>®</sup> AbSeq and Sample Tag libraries separately can help ensure each library receives enough data during sequencing due to size discrepancies. However, these libraries can be successfully sequenced in combination with WTA mRNA, TCR, and BCR libraries.
- All libraries derived from the same cartridge can be indexed with the same Illumina indices or reverse index primer from the BD Rhapsody™ reagents. The primary analysis pipeline can differentiate the library types (for example, WTA mRNA vs TCR) based on their structure and sequences. Demultiplexing statistics are reported from the pipeline, but should these statistics be desired prior to running the pipeline, then unique indices will be required for each library.

For additional support with pooling and sequencing, reach out to your local Field Application Specialist (FAS) or [scomix@bdscomix.bd.com](mailto:scomix@bdscomix.bd.com).

### Sequencing flow cell loading and PhiX concentrations

#### Quantifying libraries

Calculate the molar concentration of WTA, Sample Tag, BD<sup>®</sup> AbSeq, and TCR/BCR libraries using Qubit™ quantitation concentration (ng/μL) and average Bioanalyzer size (200–1000 bp). For TCR/BCR libraries, the expected Qubit™ concentration should be >1.5 ng/μL. Use the calculated molar concentrations to pool libraries.

#### BD<sup>®</sup> AbSeq/Sample Tag libraries

For a NextSeq High or Mid Output and MiniSeq High or Mid Output runs, load the flow cell at a concentration between 1.8–2.2 pm with 1% PhiX for a sequencing run.

**Required parameters**

Parameter	Requirement
Platform	Illumina: 150 cycle kit
Paired-end reads	Minimum of 51 × 71 paired read length
PhiX	Required (1%)
Analysis	See the <i>BD® Single-Cell Multiomics Bioinformatics Handbook</i>

**WTA and TCR/BCR libraries (with or without BD® AbSeq/Sample Tag)**

For a NextSeq High or Mid Output and MiniSeq High or Mid Output runs, load the flow cell at a concentration between 1.4–1.8 pm with 3% PhiX. For other sequencers follow Illumina recommendations for loading concentration and use 3% PhiX.

**Note:** If using less than 10,000 reads/cell for the WTA library, increase PhiX percentage to 5–10% to account for lower library diversity.

Set up sequencing run on Illumina® BaseSpace. Enter the pooled libraries as one sample if libraries were made with the same Library Forward primer but with different i7 indices.

**Required parameters**

Parameter	Requirement
Platform	Illumina: 300 cycle kit
Paired-end reads	Minimum: 65 × 150 <sup>a</sup> paired read length Recommended: 150 × 150 <sup>a</sup> paired read length
PhiX	Required (3%)
Analysis	See the <i>BD® Single-Cell Multiomics Bioinformatics Handbook</i>

a. For optimal assembly, use the sequencing configuration 85 × 215 paired read length.

**Sequencing analysis pipeline**

Contact customer support at [scomix@bdscmix.bd.com](mailto:scomix@bdscmix.bd.com) for access to the latest whole transcriptome sequencing analysis pipeline.

## Appendix

---

### Illumina Index 1 (i7) sequences

Library Reverse Primer	Sequence
1	GCTACGCT
2	CGAGGCTG
3	AAGAGGCA
4	GTAGAGGA

### Illumina Index 2 (i5) sequence

Primer	Sequence
Library Forward Primer	TATAGCCT

### Mouse T cell PCR1 primers

Primer name	Primer sequence (5'–3')
Ms_TRAC_N1	TTTTCGGCACATTGATTTGGGAG
Ms_TRBC_N1	CTCAGGCAGTAGCTATAATTGCT
Ms_TRDC_N1	CAATCTTCTGGATGATCTGAGACT
Ms_TRGC1- TRGC2_N1	GGAAAGAACTTTTCAAGGAGACAAAGG

### Mouse T cell PCR2 primers

Primer name	Primer sequence (5'–3')
Ms_TRAC_N2	AGGTTCTGGGTTCTGGATGT
Ms_TRBC_N2	CAATCTGCTTTTGATGGCTC
Ms_TRDC_N2	GTAGAAATCTTTCACCAGACAAGC
Ms_TRGC1- TRGC2_N2	TTGGGGGAAATGTCTGCA
Ms_TRGC4_N2	ATAGTAGGCTTGGGAGAAAAGTCTGA

**Mouse B cell PCR1 primers**

Primer name	Primer sequence (5'–3')
Ms_IGHA_N1	AACTGGCTGCTCATGGTGTA
Ms_IGHD_N1	AAGTGTGGTTGAGGTTCACTCTG
Ms_IGHE_N1	GAAGTTCACAGTGCTCATGTTC
Ms_IGHG1_N1	CAGAGTGTAGAGGTCAGACT
Ms_IGHG2A- IGHG2C_N1	TCGAGGTTACAGTCACTGAG
Ms_IGHG2B_N1	GATCCAGAGTTCCAAGTCACAG
Ms_IGHG3_N1	TACGTTGCAGATGACAGTCT
Ms_IGHM_N1	TGGATGACTTCAGTGTGTTCTG
Ms_IGKC_N1	TGTAGGTGCTGTCTTTGCTG
Ms_IGLC1_N1	CTGTAAGTCTATGCCTTTCCC
Ms_IGLC2-IGLC3_ N1	TTGGTGGGATTTGAAGTGTCC

**Mouse B cell PCR2 primers**

Primer name	Primer sequence (5'–3')
Ms_IGHA_N2	TGTCAGTGGGTAGATGGTGG
Ms_IGHD_N2	CTGACTTCCAATTACTAAACAGCC
Ms_IGHE_N2	TAGAGCTGAGGGTTCCTGATAG
Ms_IGHG1_N2	CAGTGGATAGACAGATGGGGGT
Ms_IGHG2A- IGHG2C_N2	ATGGGGCTGTTGTTTTGG
Ms_IGHG2B_N2	GTGGATAGACTGATGGGGGTGTT
Ms_IGHG3_N2	AGGGAAGTAGCCTTTGACAAG
Ms_IGHM_N2	GACATTTGGGAAGGACTGACTC
Ms_IGKC_N2	AGATGTAACTGCTCACTGGATG
Ms_IGLC1_N2	GTTAGTCTCGAGCTCTTCAGA
Ms_IGLC2-IGLC3_N2	CAGTGTGGCTTTGTTTTCT

## Contact Information

---

**Becton, Dickinson and Company**

**BD Biosciences**

155 North McCarthy Boulevard

Milpitas, California 95035 USA

[bdbiosciences.com](http://bdbiosciences.com)

[scomix@bdscomix.bd.com](mailto:scomix@bdscomix.bd.com)